

axiom[®]

NEW



**PRODUCT
CATALOG.**

ANTHOGYR 
A Straumann Group Brand

TABLE OF CONTENTS.

AXIOM® PROSTHETIC OVERVIEW	4
AXIOM® IMPLANTS	
Axiom® Bone Level	14
Axiom® Tissue Level	16
HEALING AND TEMPORIZATION	
Healing screws	20
Healfit® SH anatomical and scannable healing screws	22
Temporization	27
AXIOM® PROSTHESIS	
Impression at implant level	28
Labside solutions – X-Base®, CEREC®* compatible ti-bases and Pre-milled	32
Multi-Unit range	40
Esthetic and Standard Titanium abutments for cemented prosthesis	48
Removable overdentures	56
AxiIN® Simeda® CAD/CAM zirconia tooth with Angulated Access	64
Simeda® CAD/CAM plural prosthesis	68
INSTRUMENTS	
Axiom® surgery kit	70
MG kit - Dense bone kit - Countersinks	76
Drill stop kit	82
Guided surgery - INTEGRAL - FIRST DRILL	84
Axiom® prosthesis kit	92
Laboratory kits	97
Accessories for digital workflows	98
Screws and recommended torque values	102
AXIOM® 2.8	104
MINI IMPLANT SYSTEM	112
EQUIPMENT	
Torq Control®	120
Xpert Unit®	122
Mont Blanc® range	124
USEFUL INFORMATION	126

WELCOME

*CEREC® is a brand of Sirona Dental Systems GmbH

SINGLE IMPLANT SUPPORTED RESTORATIONS

TEMPORARY PROSTHESIS

AA AxIN® temporary copings
Screw-retained
with Angulated Access


BONE LEVEL
TISSUE LEVEL



→ See Page 65

Temporary abutments
Screw-retained

BONE LEVEL
TISSUE LEVEL

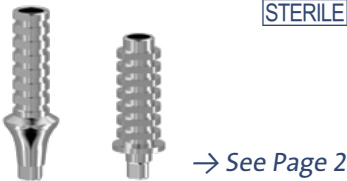

→ See Page 27

FINAL PROSTHESIS

CAD/CAM

AA AxIN® customized
screw-retained teeth
with Angulated Access

AxIN®



Simeda®

BONE LEVEL
TISSUE LEVEL

→ See Page 64

Customized screw-retained teeth
Customized abutments
for cemented restorations



Simeda®

BONE LEVEL
TISSUE LEVEL

X-Base® titanium bases
With straight or Angulated Access



LABSIDE

BONE LEVEL
TISSUE LEVEL

→ See Page 32

Pre-milled range
Customized abutments
for cemented restorations



LABSIDE

BONE LEVEL
TISSUE LEVEL

→ See Page 39

CEREC® compatible titanium bases*



CHAIRSIDE

BONE LEVEL
TISSUE LEVEL

→ See Page 38

CONVENTIONAL

Esthetic titanium abutments
for cemented restorations



→ See Page 48

Standard titanium abutments
for cemented restorations



→ See Page 52

BONE LEVEL

Components presented on these pages are not exhaustive.
To see the entire range, please refer to the page numbers indicated.

AA Angulated Access

*CEREC® is a brand of Sirona Dental Systems GmbH

PLURAL IMPLANT SUPPORTED RESTORATIONS

ABUTMENTS

TEMPORARY PROSTHESIS

FINAL PROSTHESIS

Multi-Unit abutments

BONE LEVEL

TISSUE LEVEL

→ See Page 40

Multi-Unit secondary components are the same for Axiom® BL and Axiom® TL

Multi-Unit temporary abutments with straight or Angulated Access

BONE LEVEL

TISSUE LEVEL

→ See Page 44 and 69

Temporary abutment

TISSUE LEVEL

→ See Page 27

CONVENTIONAL

Esthetic titanium abutments for cemented restorations

BONE LEVEL

→ See Page 48

Standard titanium abutments for cemented restorations

BONE LEVEL

→ See Page 52

CAD/CAM

CAD/CAM bridge, suprastructure or bar on Multi-Unit

BONE LEVEL

TISSUE LEVEL

→ See Page 32

CONVENTIONAL

Multi-Unit castable copings

BONE LEVEL

TISSUE LEVEL

→ See Page 45

X-Base® titanium base for plural, directly on Axiom® TL With straight or Angulated Access

TISSUE LEVEL

→ See Page 32

Pre-milled range for customized abutments

BONE LEVEL

TISSUE LEVEL

→ See Page 39

SOLUTIONS FOR REMOVABLE OVERDENTURES

Novaloc® - Locator®

BONE LEVEL

TISSUE LEVEL

→ See Page 56

axiom[®]

DRIVEN BY
ANTHOGYR **EXPERTISE.**



80 years of French expertise,
shaping the future of implant dentistry.

- Precision engineered in the French Alps
- Innovation driven by simplicity, born to seamlessly integrate into digital workflows
- Together, beyond products.

axiom[®]


THE POWER OF
SIMPLICITY.

High performance
designs to take up your
clinical challenges.

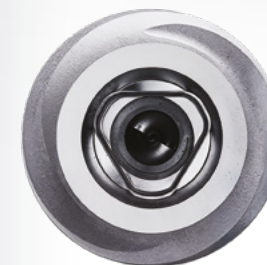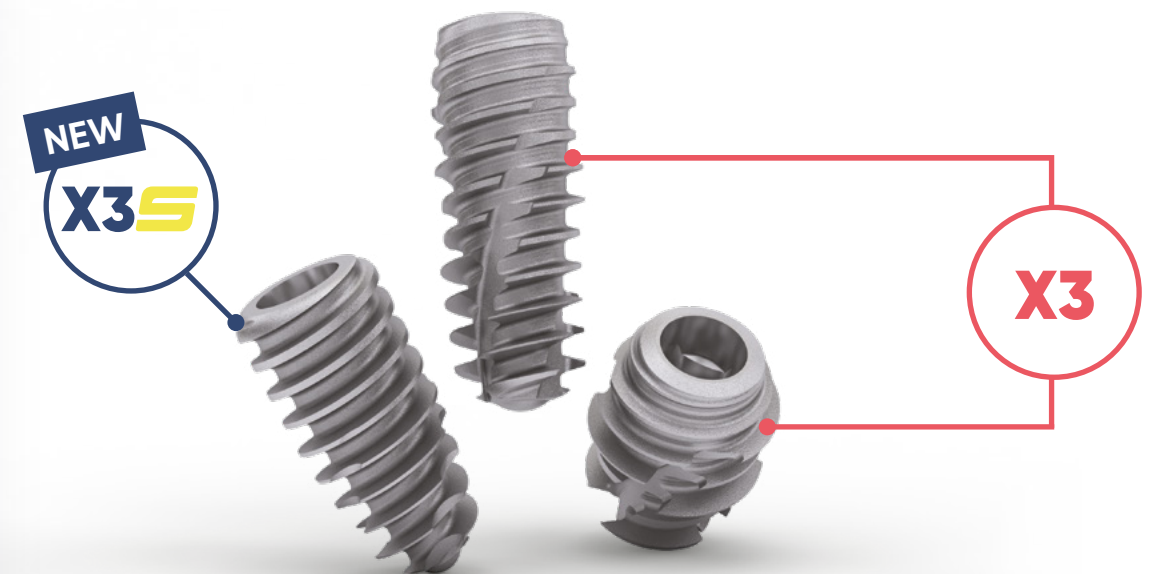

1 UNIQUE
prosthetic
platform

1 single COMPACT
surgical kit

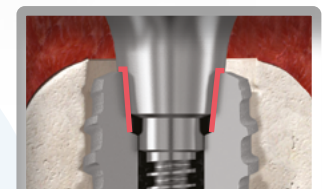

1 SEALED conical
connection for long
term bone stability⁽¹⁾

Seamlessly **INTEGRATED** into the digital workflow.

(1) Esposito et al.. Ten-year post loading outcomes of subcrestally placed implant with internal conical connections at 0.5 mm versus 1.5 mm depth: a multicentre within-person randomised controlled trial. Clinical Oral Implants research. 2026; 0:1-17.

HIGH PERFORMANCE IMPLANT DESIGNS FOR SMART BONE MANAGEMENT

NEW

axiom®X3S

Enhanced anchorage in demanding clinical situations⁽¹⁾

STRAIGHT NECK

Cortical anchorage

WIDE THREADS

Lateral anchorage

TAPERED DESIGN

Bone compression

3 BI-ACTION CUTTING FLUTES

Guided insertion

1. ACTIVE CUTTING EDGE


Apical anchorage

2. BONE PRESERVATION AND COLLECTION

For progressive stabilization

SLIM AND ROUNDED APEX

Engagement in narrow osteotomies
Designed to protect sinus membrane

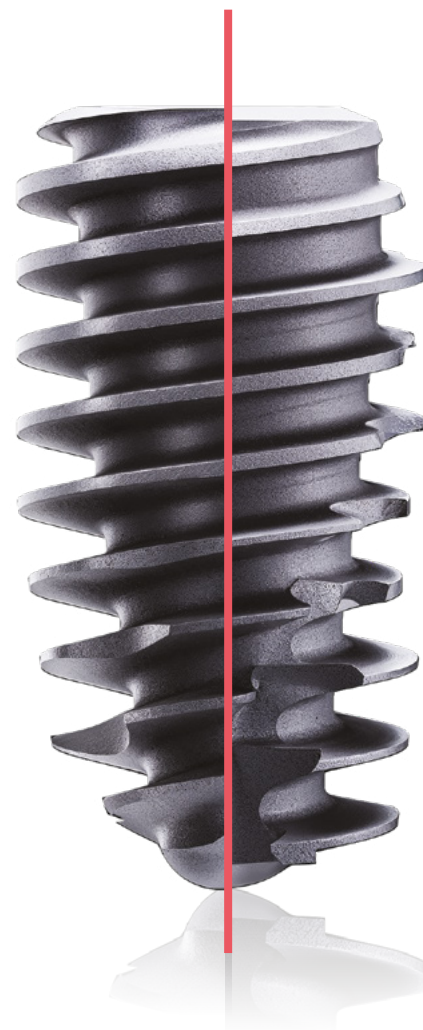

axiom®X3

Innovative implant design for universal bone anchorage and bone preservation

REDUCED NECK DIAMETER

Less cortical stress

ALTERNATIVE THREAD LENGTH

Bone preservation
Lower insertion efforts
All bone densities

SLIM THREADS

Easier insertion

TAPERED PROFILE

Optimal primary stability
in soft bone

3 BI-DIRECTIONAL CUTTING FLUTES

Dynamic engagement all over the implant height
Insertion axis correction

ROUNDED AND THIN APEX

Sinus membrane protection
Efficient engagement in post extractional situations

SINCE 2009



axiom®REG

Parallel wall, versatile implant with progressive anchorage

STRAIGHT THREADED NECK

Designed for optimized primary bone anchorage

CYLINDRICAL - CONICAL SHAPE

Progressive compression
Controlled and guided insertion

ASYMMETRICAL THREADS

Self-tapping
Progressive anchorage

ROUNDED APEX

Gentle shape when approaching the sinus



axiom®PX

Apically tapered implant, for immediate post-extraction and low-density bone

REVERSED CONICAL NECK*

Preservation of cortical bone volume

TAPERED DESIGN

Designed for initial stability in low-bone density

DEEP SYMMETRICAL DOUBLE THREAD

Guiding effect, fast bone penetration
Self-drilling and self-tapping

BITING APEX

Designed for optimal anchorage on the apical portion

(1) Axiom® X3S reached higher primary stability than Axiom® X3 in D3 - D4 bone during International Product Launch Survey. Data on file.

CHOOSE YOUR
CLINICAL APPROACH

BONE LEVEL

Endosseous implant designed surgical and prosthetic flexibility.



SINGLE DIAMETER PROSTHETIC CONNECTION

Prosthetic convenience



SEALED CONICAL CONNECTION

Designed for biological integration.



PLATFORM SWITCHING

Soft tissue stabilization.



INTUITIVE TRILOBE

Guided, accurate abutment placement.

CHOICE OF DESIGNS X3 X3 PX REG

TISSUE LEVEL

Transgingival implant designed for biological process respect.




TISSUE FAVORED DESIGN

Anatomical neck profile for soft tissue maximization and preservation.



6 EMERGENCE PROFILES

2 platforms

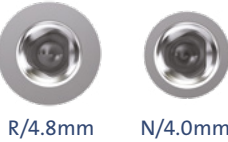

3 gingival heights



A TRUSTED REFERENCE FOR COMPROMISED PATIENTS

High performance in real-world conditions. 98.9% implant survival rate⁽¹⁾

CHOICE OF DESIGNS X3 PX REG

(1) Chacun D, Lefebvre L, Langlois-Jacques C, Dard M, Gritsch K, Grosgeat B. A Multicenter Cohort Study on 301 Tissue-Level Implants: Cumulative Implant Survival Rate and Marginal Bone Level Change up to 4.5 Years. Int J Oral Maxillofac Implants. 2024 Apr 24;39(2):224-234. doi: 10.11607/jomi.10141. PMID: 38265484.

Axiom® Bone Level implants

axiom®X3

STERILE

DIAMETER

Ø 3.4 mm ●

Ø 4.0 mm ●

Ø 4.6 mm ○

Ø 5.2 mm ●

Ø 5.8mm ●

Ø 6.4 mm ●

COVER SCREW INCLUDED

OPIM100*

LENGTH

6.5 mm

8 mm

10 mm

12 mm

14 mm

16 mm

18 mm

6.5 mm	-	XT40060	XT46060	XT52060	XT58060	XT64060
8 mm	XT34080	XT40080	XT46080	XT52080	XT58080	XT64080
10 mm	XT34100	XT40100	XT46100	XT52100	XT58100	XT64100
12 mm	XT34120	XT40120	XT46120	XT52120	XT58120	XT64120
14 mm	XT34140	XT40140	XT46140	-	-	-
16 mm	XT34160	XT40160	-	-	-	-
18 mm	XT34180	XT40180	-	-	-	-

NEW

STERILE

DIAMETER

Ø 3.4 mm ●

Ø 4.0 mm ●

Ø 4.6 mm ○

Ø 5.2 mm ●

COVER SCREW INCLUDED

OPIM100*

LENGTH

6.5 mm

8 mm

10 mm

12 mm

14 mm

6.5 mm	-	AT40060	AT46060	AT52060
8 mm	AT34080	AT40080	AT46080	AT52080
10 mm	AT34100	AT40100	AT46100	AT52100
12 mm	AT34120	AT40120	AT46120	AT52120
14 mm	AT34140	AT40140	AT46140	-

All implant visuals are shown in 10 mm length.

TOP REFERENCES ARE SHOWN IN CORAL

BONE HEIGHT	0.5 mm	1.0 mm	1.5 mm	2.0 mm	2.5 mm
Cover screws					
	OPIM100*	OPIM110	OPIM115	OPIM120	OPIM125

*Manual surgical wrench

OPCS100

Osstell ISQ SmartPeg Type to be used for Axiom® BL: Type 38

*included with each Axiom® Bone Level implant



Axiom® Tissue Level implants

STERILE		axiom [®] REG				axiom [®] PX			
DIAMETER		Ø 3.4 mm ●	Ø 4.0 mm ●	Ø 4.6 mm ○	Ø 5.2 mm ●	Ø 3.4 mm ●	Ø 4.0 mm ●	Ø 4.6 mm ○	Ø 5.2 mm ●
Manual ≤ 10 N.cm	Cover screw included								
	OPIM100*								
LENGTH									
6.5 mm		-	OP40060	OP46060	OP52060	-	PX40060	PX46060	PX52060
8 mm		OP34080	OP40080	OP46080	OP52080	PX34080	PX40080	PX46080	PX52080
10 mm		OP34100	OP40100	OP46100	OP52100	PX34100	PX40100	PX46100	PX52100
12 mm		OP34120	OP40120	OP46120	OP52120	PX34120	PX40120	PX46120	PX52120
14 mm		OP34140	OP40140	OP46140	OP52140	PX34140	PX40140	PX46140	-
16 mm		OP34160	OP40160	-	-	PX34160	PX40160	-	-
18 mm		OP34180	OP40180	-	-	PX34180	PX40180	-	-

All implant visuals are shown in 10 mm length.

*Manual surgical wrench	
OPCS100	

Osstell ISQ SmartPeg Type to be used for Axiom® TL: Type 66

Axiom® Tissue Level implants

		axiomX3®						PLATFORM N/4.0 mm	
		PLATFORM R/4.8 mm							
		Ø 3.4 mm	Ø 4.0 mm	Ø 4.6 mm	Ø 5.2 mm	Ø 5.8 mm	Ø 6.4 mm	Ø 3.4 mm	Ø 4.0 mm
Manual ≤ 10 N.cm	Sealing screw included								
	Manual surgical wrench OPCS100								
IMPLANT LENGHT									
NECK HEIGHT									
6.5 mm	1.5 mm	-	TXT40060R1	TXT46060R1	TXT52060R1	TXT58060R1	TXT64060R1	-	TXT40060N1
	2.5 mm	-	TXT40060R2	TXT46060R2	TXT52060R2	TXT58060R2	TXT64060R2	-	TXT40060N2
	3.5 mm	-	TXT40060R3	TXT46060R3	-	-	-	-	TXT40060N3
8 mm	1.5 mm	TXT34080R1	TXT40080R1	TXT46080R1	TXT52080R1	TXT58080R1	TXT64080R1	TXT34080N1	TXT40080N1
	2.5 mm	TXT34080R2	TXT40080R2	TXT46080R2	TXT52080R2	TXT58080R2	TXT64080R2	TXT34080N2	TXT40080N2
	3.5 mm	TXT34080R3	TXT40080R3	TXT46080R3	-	-	-	TXT34080N3	TXT40080N3
10 mm	1.5 mm	TXT34100R1	TXT40100R1	TXT46100R1	TXT52100R1	TXT58100R1	TXT64100R1	TXT34100N1	TXT40100N1
	2.5 mm	TXT34100R2	TXT40100R2	TXT46100R2	TXT52100R2	TXT58100R2	TXT64100R2	TXT34100N2	TXT40100N2
	3.5 mm	TXT34100R3	TXT40100R3	TXT46100R3	-	-	-	TXT34100N3	TXT40100N3
12 mm	1.5 mm	TXT34120R1	TXT40120R1	TXT46120R1	TXT52120R1	TXT58120R1	TXT64120R1	TXT34120N1	TXT40120N1
	2.5 mm	TXT34120R2	TXT40120R2	TXT46120R2	TXT52120R2	TXT58120R2	TXT64120R2	TXT34120N2	TXT40120N2
	3.5 mm	TXT34120R3	TXT40120R3	TXT46120R3	-	-	-	TXT34120N3	TXT40120N3
14 mm	1.5 mm	TXT34140R1	TXT40140R1	TXT46140R1	-	-	-	TXT34140N1	TXT40140N1
	2.5 mm	TXT34140R2	TXT40140R2	TXT46140R2	-	-	-	TXT34140N2	TXT40140N2
	3.5 mm	TXT34140R3	TXT40140R3	TXT46140R3	-	-	-	TXT34140N3	TXT40140N3
16 mm	1.5 mm	TXT34160R1	TXT40160R1	-	-	-	-	TXT34160N1	TXT40160N1
	2.5 mm	TXT34160R2	TXT40160R2	-	-	-	-	TXT34160N2	TXT40160N2
	3.5 mm	TXT34160R3	TXT40160R3	-	-	-	-	TXT34160N3	TXT40160N3

All implant visuals are shown in 10 mm length, and 2.5 mm neck height.
Osstell ISQ SmartPeg Type to be used for Axiom® TL: Type 66

Calibrating film / X-Ray template

Axiom X3® Tissue Level calibrating film: TXTFC_NOT

Axiom® Tissue Level implants

Manual
10 N.cm

Sealing screw included

TSS-R0

TSS-N0

Manual surgical wrench
OPCS100

STERILE

IMPLANT
LENGHT

NECK
HEIGHT

axiom®REG

		PLATFORM R/4.8 mm				PLATFORM N/4.0 mm	
		Ø 3.4 mm ●	Ø 4.0 mm ●	Ø 4.6 mm ○	Ø 5.2 mm ●	Ø 3.4 mm ●	Ø 4.0 mm ●
6.5 mm	1.5 mm	-	TOP40060R1	TOP46060R1	TOP52060R1	-	TOP40060N1
	2.5 mm	-	TOP40060R2	TOP46060R2	TOP52060R2	-	TOP40060N2
	3.5 mm	-	TOP40060R3	TOP46060R3	-	-	TOP40060N3
8 mm	1.5 mm	TOP34080R1	TOP40080R1	TOP46080R1	TOP52080R1	TOP34080N1	TOP40080N1
	2.5 mm	TOP34080R2	TOP40080R2	TOP46080R2	TOP52080R2	TOP34080N2	TOP40080N2
	3.5 mm	TOP34080R3	TOP40080R3	TOP46080R3	-	TOP34080N3	TOP40080N3
10 mm	1.5 mm	TOP34100R1	TOP40100R1	TOP46100R1	TOP52100R1	TOP34100N1	TOP40100N1
	2.5 mm	TOP34100R2	TOP40100R2	TOP46100R2	TOP52100R2	TOP34100N2	TOP40100N2
	3.5 mm	TOP34100R3	TOP40100R3	TOP46100R3	-	TOP34100N3	TOP40100N3
12 mm	1.5 mm	TOP34120R1	TOP40120R1	TOP46120R1	TOP52120R1	TOP34120N1	TOP40120N1
	2.5 mm	TOP34120R2	TOP40120R2	TOP46120R2	TOP52120R2	TOP34120N2	TOP40120N2
	3.5 mm	TOP34120R3	TOP40120R3	TOP46120R3	-	TOP34120N3	TOP40120N3
14 mm	1.5 mm	TOP34140R1	TOP40140R1	TOP46140R1	-	TOP34140N1	TOP40140N1
	2.5 mm	TOP34140R2	TOP40140R2	TOP46140R2	-	TOP34140N2	TOP40140N2
	3.5 mm	TOP34140R3	TOP40140R3	TOP46140R3	-	TOP34140N3	TOP40140N3
16 mm	1.5 mm	TOP34160R1	TOP40160R1	-	-	TOP34160N1	TOP40160N1
	2.5 mm	TOP34160R2	TOP40160R2	-	-	TOP34160N2	TOP40160N2
	3.5 mm	TOP34160R3	TOP40160R3	-	-	TOP34160N3	TOP40160N3

Osstell ISQ SmartPeg Type to be used for Axiom® TL: Type 66

Axiom® Tissue Level implants

Manual
10 N.cm

Sealing screw included

TSS-R0

TSS-N0

Manual surgical wrench
OPCS100

STERILE

IMPLANT
LENGHT

NECK
HEIGHT

axiom®PX

		PLATFORM R/4.8 mm				PLATFORM N/4.0 mm	
		Ø 3.4 mm ●	Ø 4.0 mm ●	Ø 4.6 mm ○	Ø 5.2 mm ●	Ø 3.4 mm ●	Ø 4.0 mm ●
6.5 mm	1.5 mm	-	TPX40060R1	TPX46060R1	TPX52060R1	-	TPX40060N1
	2.5 mm	-	TPX40060R2	TPX46060R2	TPX52060R2	-	TPX40060N2
	3.5 mm	-	TPX40060R3	TPX46060R3	-	-	TPX40060N3
8 mm	1.5 mm	TPX34080R1	TPX40080R1	TPX46080R1	TPX52080R1	TPX34080N1	TPX40080N1
	2.5 mm	TPX34080R2	TPX40080R2	TPX46080R2	TPX52080R2	TPX34080N2	TPX40080N2
	3.5 mm	TPX34080R3	TPX40080R3	TPX46080R3	-	TPX34080N3	TPX40080N3
10 mm	1.5 mm	TPX34100R1	TPX40100R1	TPX46100R1	TPX52100R1	TPX34100N1	TPX40100N1
	2.5 mm	TPX34100R2	TPX40100R2	TPX46100R2	TPX52100R2	TPX34100N2	TPX40100N2
	3.5 mm	TPX34100R3	TPX40100R3	TPX46100R3	-	TPX34100N3	TPX40100N3
12 mm	1.5 mm	TPX34120R1	TPX40120R1	TPX46120R1	TPX52120R1	TPX34120N1	TPX40120N1
	2.5 mm	TPX34120R2	TPX40120R2	TPX46120R2	TPX52120R2	TPX34120N2	TPX40120N2
	3.5 mm	TPX34120R3	TPX40120R3	TPX46120R3	-	TPX34120N3	TPX40120N3
14 mm	1.5 mm	TPX34140R1	TPX40140R1	TPX46140R1	-	TPX34140N1	TPX40140N1
	2.5 mm	TPX34140R2	TPX40140R2	TPX46140R2	-	TPX34140N2	TPX40140N2
	3.5 mm	TPX34140R3	TPX40140R3	TPX46140R3	-	TPX34140N3	TPX40140N3
16 mm	1.5 mm	TPX34160R1	TPX40160R1	-	-	TPX34160N1	TPX40160N1
	2.5 mm	TPX34160R2	TPX40160R2	-	-	TPX34160N2	TPX40160N2
	3.5 mm	TPX34160R3	TPX40160R3	-	-	TPX34160N3	TPX40160N3

TOP REFERENCES ARE SHOWN IN CORAL

All implant visuals are shown in 10 mm length.
Osstell ISQ SmartPeg Type to be used for Axiom® TL: Type 66

Sealing screws

Platform identification
by laser marking "TL N" and "TL R"

TISSUE LEVEL	
PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
 TSS-R0	 TSS-N0

18

ANTHOGRY

19

Healing screws Axiom® BL


Check our tips to choose the most adapted healing screw page 146



STERILE GINGIVAL HEIGHT	BONE LEVEL							
	Ø 3.4 mm		Ø 4.0 mm		Ø 5.0 mm		Ø 6.0 mm	
	CORONAL HEIGHT		CORONAL HEIGHT		CORONAL HEIGHT		CORONAL HEIGHT	
	2 mm	1 mm	2 mm	1 mm	2 mm	1 mm	2 mm	1 mm
0.75 mm	-	-					-	-
			OPHS400	OPHSF400	OPHS500	OPHSF500		
1.5 mm								
	OPHS310	OPHSF310	OPHS410	OPHSF410	OPHS510	OPHSF510	OPHS610	OPHSF610
2.5 mm								
	OPHS320	OPHSF320	OPHS420	OPHSF420	OPHS520	OPHSF520	OPHS620	OPHSF620
3.5 mm								
	OPHS330	OPHSF330	OPHS430	OPHSF430	OPHS530	OPHSF530	OPHS630	OPHSF630
4.5 mm*								
	OPHS340*	OPHSF340*	OPHS440*	OPHSF440*	OPHS540*	OPHSF540*	OPHS640*	OPHSF640*

LASER MARKING ON SCREW HEAD

GH (mm)	0.75	1.5	2.5	3.5	4.5
Ø 3.4					
Ø 4.0					
Ø 5.0					
Ø 6.0					

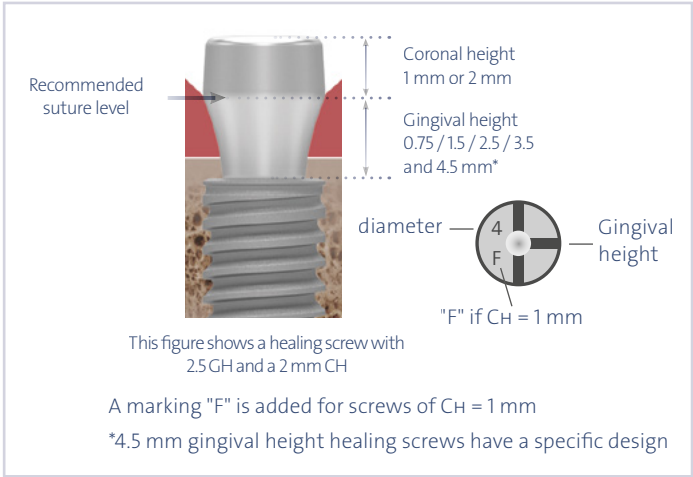

Healing screws Axiom® TL



Healing screws

STERILE Manual surgical wrench OPCS100	TISSUE LEVEL			
	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
	2 mm	4 mm	2 mm	4 mm
Conical profile				
	THS-R210	THS-R410	THS-N210	THS-N410
Straight profile		-		-
	THS-R200		THS-N200	

healfit® SH

Anatomical and scannable healing solution with biomimetic design



- 3-in-1 solution: Heal – Shape – Take the impression
- For single restorations on Axiom® Bone Level and Tissue Level
- 5 anatomical shapes based on the study of hundreds of tooth profiles
- Choice of gingival and coronal heights
- Open software solution available in major planning softwares, with a direct access to CAD libraries on Anthogyr website

SCANBODY FUNCTION

- Fully integrated digital workflow
- Serves as transfer

ZrN* COATING

- Natural appearance
- Fast and efficient scanning
- Biocompatible, with antibacterial properties

AXIOM® CONNECTION

- For Axiom® BL and TL

VISUAL IDENTIFICATION


- Visual decoding of relevant information on the occlusal surface: gingival and coronary height, implant type, trilobe orientation

CONSTANT EMERGENCE PROFILE

- Consistent with X-Base® profiles

GOOD TO KNOW


- For proper impression taking, coronal part must exceed gum line at least +1.5 mm

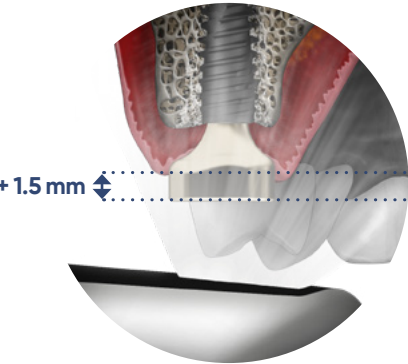

*Zirconium Nitride


5 unique shapes*

Inspired from the study of hundreds of dental profiles, to offer an anatomical fit across all sectors

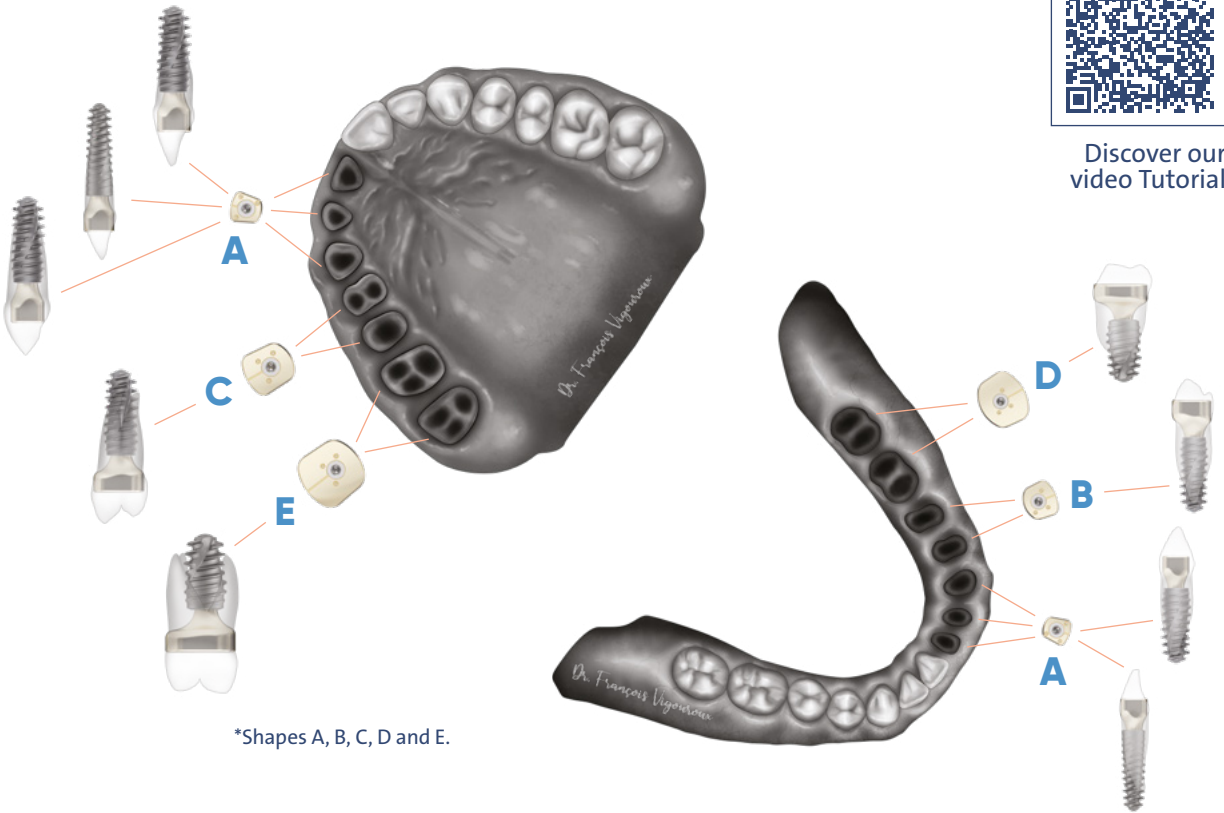

*Shapes A, B, C, D and E.

A choice of gingival and coronal heights

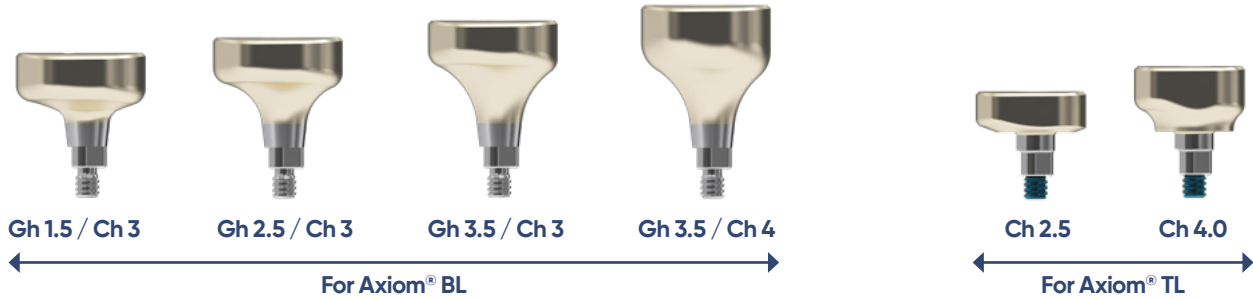

Dedicated instrument


Height gauge for HealFit® SH BL and TL



OPJCSHS

*Gh: Gingival height - *Ch: Coronary height


Anatomical and scannable healing screws

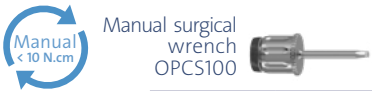

BONE LEVEL				TISSUE LEVEL			
SHAPE	Gh*	Ch*	FOR SINGLE	SHAPE	Ch*	PLATFORM R/4.8 MM	PLATFORM N/4.0 MM
A	1.5 mm	3 mm	OPSHSA13 	A	2.5 mm	-	TSHSA-N200 
	2.5 mm	3 mm	OPSHSA23 				
	3.5 mm	3 mm	OPSHSA33 		4.0 mm	-	TSHSA-N400 
		4 mm	OPSHSA34 				
B	1.5 mm	3 mm	OPSHSB13 	B	2.5 mm	TSHSB-R200 	TSHSB-N200 
	2.5 mm	3 mm	OPSHSB23 				
	3.5 mm	3 mm	OPSHSB33 		4.0 mm	TSHSB-R400 	TSHSB-N400 
		4 mm	OPSHSB34 				
C	1.5 mm	3 mm	OPSHSC13 	C	2.5 mm	TSHSC-R200 	TSHSC-N200 
	2.5 mm	3 mm	OPSHSC23 				
	3.5 mm	3 mm	OPSHSC33 		4.0 mm	TSHSC-R400 	TSHSC-N400 
		4 mm	OPSHSC34 				
D	1.5 mm	3 mm	OPSHSD13 	D	2.5 mm	TSHSD-R200 	-
	2.5 mm	3 mm	OPSHSD23 				
	3.5 mm	3 mm	OPSHSD33 		4.0 mm	TSHSD-R400 	-
		4 mm	OPSHSD34 				
E	1.5 mm	3 mm	OPSHSE13 	E	2.5 mm	TSHSE-R200 	-
	2.5 mm	3 mm	OPSHSE23 				
	3.5 mm	3 mm	OPSHSE33 		4.0 mm	TSHSE-R400 	-
		4 mm	OPSHSE34 				

*Gh: Gingival height - *Ch: Coronary height

Healfit® SH for Axiom® BL

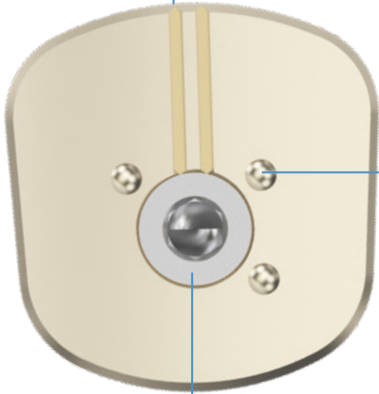

Coronary height indicator

- 1 groove = Ch 3 mm
- 2 grooves = Ch 4 mm

💡 Always oriented on vestibular side!

Gingival height indicator

- 1 dot = Gh 1.5mm
- 2 dots = Gh 2.5 mm
- 3 dots = Gh 3.5 mm


Grey screw = Axiom® BL

GOOD TO KNOW !

Healfit® SH body and screw delivered unassembled



Healfit® SH for Axiom® TL

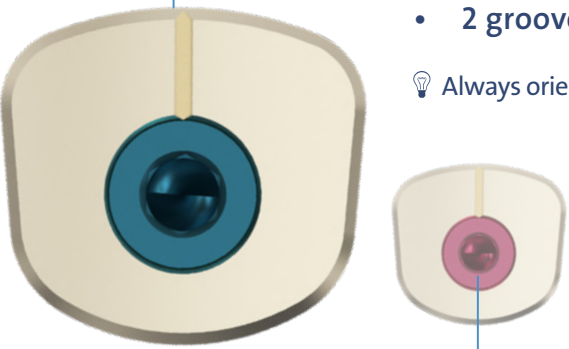

Coronary height indicator

- 1 groove = Ch 2.5 mm
- 2 grooves = Ch 4 mm


💡 Always oriented on vestibular side!

Blue screw = Axiom® TL R platform
Pink screw = Axiom® TL N platform

GOOD TO KNOW !

Healfit® SH body and screw delivered assembled



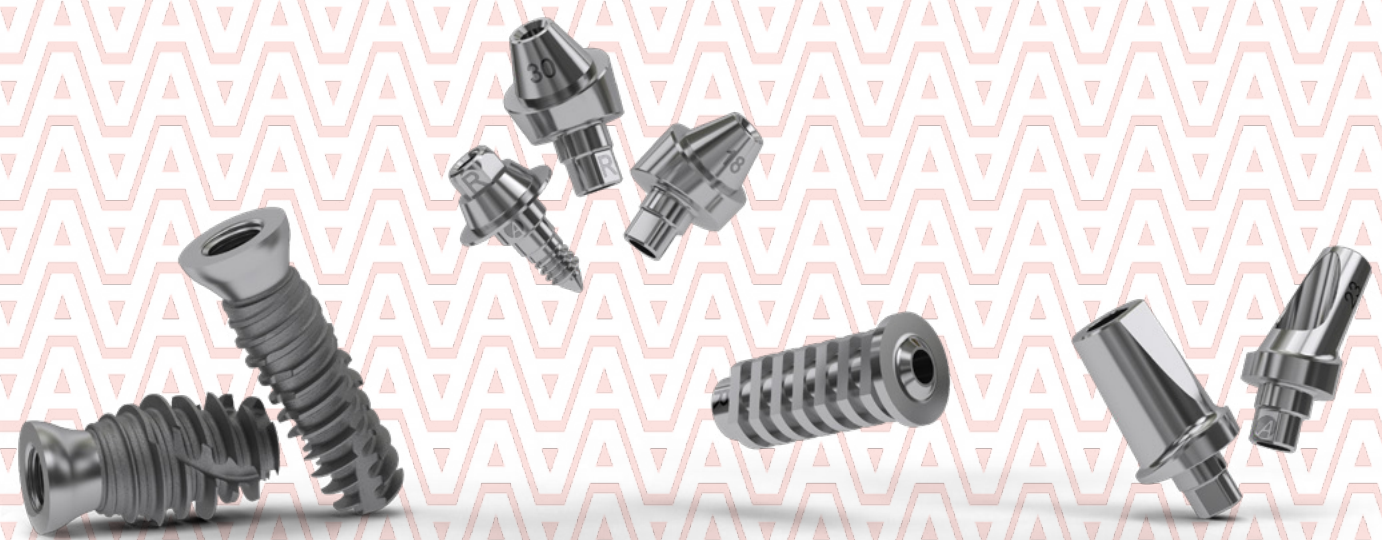

TEMPORARY ABUTMENTS

Temporary restoration

- Retention grooves
- Reworkable

CH = 9 mm

Ø

GH

STERILE

Axiom® BL

Supplied with

25 N.cm

Manual

M1.6 titanium screw OPTS161

Laboratory screw OPTS162

Platform identification by laser marking «N» and «R»

R

Axiom® TL

Supplied with

25 N.cm

Manual

TS161

Laboratory screw TS163 - TS162

25 N.cm

Axiom® TL

Supplied with

25 N.cm

Manual

TS161P

Laboratory screw TS162P-2 - TS163P-2

For Bone Level, removable with the gripper INEXPS/L M1.6 titanium screw included OPTS161

GINGIVAL HEIGHT

 For Bone Level, removable with the gripper INEXPS/L M1.6 titanium screw included OPTS161 GINGIVAL HEIGHT	BONE LEVEL				TISSUE LEVEL			
	Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
					Single (indexed)	Plural (not indexed)	Single (indexed)	Plural (not indexed)
								
H 0.75 mm	-	OPTP400	OPTP500	-				
H 1.5 mm	OPTP310	OPTP410	OPTP510	OPTP610				
H 2.5 mm	OPTP320	OPTP420	OPTP520	OPTP620				
H 3.5 mm	OPTP330	OPTP430	OPTP530	OPTP630				
H 4.5 mm	OPTP340	OPTP440	OPTP540	OPTP640	TC100-R	TC100-R-P	TC100-N	TC100-N-P

IMPRESSION RANGE

A complete range to cover all technics

- Axiom® BL - One transfer - One analog for all implant diameters
- Impression rings for pick-up transfers on Axiom® BL
- Axiom® TL - Color coding of platforms R and N
- Digital or conventional impression

● PLATFORM R/4.8 mm
● PLATFORM N/4.0 mm



Components for conventional impression taking

BONE LEVEL			TISSUE LEVEL			
Manual Pop-in transfer Screw included	SHORT	OPPI100S (X1) OPPI100S-4 (X4)	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
			Single (indexed)	Plural (not indexed)*	Single (indexed)	Plural (not indexed)*
Manual Pick-up-transfer 2 screws included	SHORT	OPPI100 (X1) OPPI100-4 (X4)	TT200-R	ILT200-R-4 (X4)	TT200-N	ILT200-N-4 (X4)
	LONG	OPPI100L (X1) OPPI100L-4 (X4)	TT200L-R	ILT100-R ILT100-R-4 (X4)	TT200L-N	ILT100-N ILT100-N-4 (X4)
Hexagonal screwdriver	COMPATIBLE GINGIVAL HEIGHT	Ø 4 mm	Ø 5 mm			
	1.5 mm	OPROFIL410 (X4)	OPROFIL510 (X4)			
	2.5 mm / 4.5 mm	OPROFIL420 (X4)	OPROFIL520 (X4)			
	3.5 mm	OPROFIL430 (X4)	OPROFIL530 (X4)			
Implant analog	Laboratory screw included OPIA100 (X1) - OPIA100-4 (X4)		TA100-R (X1) TA100-R-4 (X4)		TA100-N (X1) TA100-N-4 (X4)	

Hexagonal screwdriver.

Components for digital impression from an Intra-Oral Scanner

Manual
5N.cm

	BONE LEVEL	TISSUE LEVEL			
		PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
		SINGLE	PLURAL*	SINGLE	PLURAL*
Digital transfer Screw included	 152-27-DT 12.5 mm	 156-02-DT	 156-02-DT-IL	 156-01-DT	 156-01-DT-IL
Metal digital transfer Screw included	 152-27-MDT 13.5 mm	 156-02-S-MDT		 156-01-S-MDT	
Analog for printed model	 Laboratory screw included 152-27-PA 152-27-PA-4 (X4)	 156-02-PA 156-02-PA-4 (X4)		 156-01-PA 156-01-PA-4 (X4)	

Hexagonal screwdriver.

Scan-Adapters for digitization of conventional impression

Manual
5N.cm

	BONE LEVEL	TISSUE LEVEL	
		SINGLE R AND N PLATFORMS	PLURAL* R AND N PLATFORMS
Scan-Adapter Laboratory screw included	 152-27-SAA	 156-0X-SAA	 156-0X-SAO
Scan-Adapter screwdriver	 SATOOL-01		



X-BASE® RANGE

Ti-base range for customized screw-retained restorations, produced in dental laboratory

- For Axiom® Bone Level, Tissue Level and Multi-Unit
- For single and plural restorations
- Straight and Angulated Access, up to 25°
- Choice of gingival heights, coronary heights and platforms
- LaserGrip® surface treatment allowing time saving and secure bonding
- Freedom of CAD, CAM and milling solution
- Delivered with the associated definitive screw



Discover our video Tutorials

X-Base® with straight access

Ti-bases for Axiom® BL are extractable with INEXPS/L or INEXPHS/L

BONE LEVEL					TISSUE LEVEL					
Gh*	Ch*	FOR SINGLE			Ch*	FOR SINGLE		FOR PLURAL		
		Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm		PLATFORM		Ch*	PLATFORM	
						R/4.8 mm	N/4.0 mm		R/4.8 mm	N/4.0 mm
1.5	4				4			4		
	6									
2.5	4				6			5		
	6									
3.5	4				6			5		
	6									
Cast coping					6			5		
Definitive screw										
Short laboratory screw										
Long laboratory screw										


Hexagonal screwdriver.

X-Base® with Angulated Access for SINGLE

Ti-bases for Axiom BL are extractible with INEXPHS/L



BONE LEVEL					TISSUE LEVEL	
Gh*	Ch*	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	Ch*	FOR SINGLE
						PLATFORM
						R/4.8 mm N/4.0 mm
1.5	4	 OPFLEX414-AA / OPFLEX414-AAU	 OPFLEX514-AA OPFLEX514-AAU	 OPFLEX614-AA OPFLEX614-AAU	4	 TFLEX-R4-SAA TFLEX-R4-SAAU
	6	 OPFLEX416-AA OPFLEX416-AAU	 OPFLEX516-AA OPFLEX516-AAU	 OPFLEX616-AA OPFLEX616-AAU		
2.5	4	 OPFLEX424-AA OPFLEX424-AAU	 OPFLEX524-AA OPFLEX524-AAU	 OPFLEX624-AA OPFLEX624-AAU	4	 TFLEX-R4-SAA TFLEX-R4-SAAU
	6	 OPFLEX426-AA OPFLEX426-AAU	 OPFLEX526-AA OPFLEX526-AAU	 OPFLEX626-AA OPFLEX626-AAU		
3.5	4	 OPFLEX434-AA OPFLEX434-AAU	 OPFLEX534-AA OPFLEX534-AAU	 OPFLEX634-AA OPFLEX634-AAU	6	 TFLEX-R6-SAA TFLEX-R6-SAAU
	6	 OPFLEX436-AA OPFLEX436-AAU	 OPFLEX536-AA OPFLEX536-AAU	 OPFLEX636-AA OPFLEX636-AAU		
Cast coping AA		 OPFLEXC4-AA25	 OPFLEXC5-AA25	 OPFLEXC6-AA25	6 AA 25°	 TFLEXC-R-AA25 TFLEXC-N-AA25

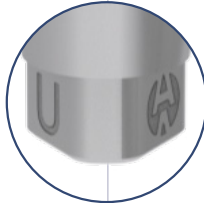

«U» version = opposite side trilobe


X-Base® with Angulated Access for PLURAL

Ti-bases for Axiom BL are extractible with INEXPHS/L



BONE LEVEL					TISSUE LEVEL	
Gh	Ch	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	Ch*	PLATFORM
						R/4.8 mm N/4.0 mm
1.5	5	 OPFLEX415-PAA	 OPFLEX515-PAA	 OPFLEX615-PAA	4	 TFLEX-R4-PAA TFLEX-N4-PAA
	5	 OPFLEX425-PAA	 OPFLEX525-PAA	 OPFLEX625-PAA		
3.5	5	 OPFLEX435-PAA	 OPFLEX535-PAA	 OPFLEX635-PAA	4 AA 25°	 TFLEXCRP-AA25 TFLEXCNP-AA25

X-Base for bridge on Axiom® BL

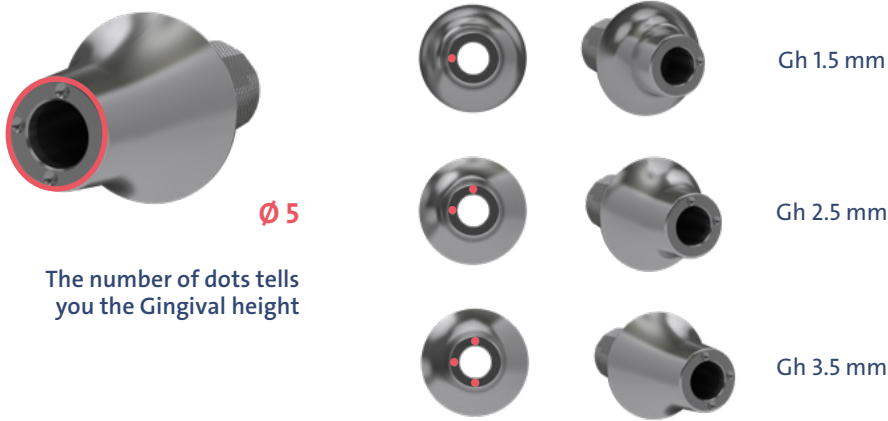

For implant-supported bridges on Axiom® BL, we recommend the use of Multi-Unit abutments combined with X-Base, offering enhanced restorative versatility.

Screw for X-Base® with Angulated Access

	Gh*	BONE LEVEL	TISSUE LEVEL	
Definitive screw AA	1.5		SINGLE  TFLEXS-AA	PLURAL  TFLEXS-PAA
		OPFLEXS1-AA		
	2.5			
	OPFLEXS2-AA			
3.5				
	OPFLEXS3-AA			
Short laboratory screw AA		BONE LEVEL	TISSUE LEVEL	
	1.5		SINGLE  TFLEXSL-AA	PLURAL  TFLEXSL-PAA-2
	2.5			
	OPFLEXSL2-AA			
3.5				
OPFLEXSL3-AA				

Instrument for Angulated Access

Angulated Access screw gripper	 AATOOL AATOOL-4 (X4)		 AATOOL-M
Ball mandrel	 INBM100S-3 (X3) 20.5 mm	 INBM100L-3 (X3) 27 mm	 INBM100XL-3 (X3)* 35 mm
Ball wrench	 INBW100S-3 (X3) 23 mm	 INBW100L-3 (X3) 29 mm	 INBW100XL-3 (X3)* 37 mm

X-Base® for plural screw-retained restorations on Multi-Unit abutments

Definitive screw included 	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
X-Base® Multi-Unit - Straight	Coronal heigh: 5 mm  MUFLEX-5	Coronal heigh: 4 mm  MUNFLEX-4		
Cast coping	 MUFLEXC		 MUNFLEXC	
Definitive screw	 MU140			
Laboratory screw	Short  MU141	Long  MUT101	Longer  MUT102	
Definitive screw included 	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
	Coronal heigh: 5 mm  MUFLEX-5-AA	Coronal heigh: 4 mm  MUNFLEX-4-AA		
Cast coping AA 25°	 MUFLEXC-AA25		 MUNFLEXC-AA25	
Definitive screw	 MUAA140			
Short laboratory screw	 MUAA142-4 (X4)			

CEREC®* COMPATIBLE TI-BASE RANGE

For customized single-unit screw-retained restorations, produced in dental office

A consistent range for Axiom® BL and TL

- Axiom® BL: 3 gingival heights
- Axiom® TL: 2 platforms
- Compatible with Dentsply Sirona® CEREC® digitalization and milling solutions



For Bone Level, removable with the gripper INEXPS/L or INEXPHS/L

GINGIVAL HEIGHT	BONE LEVEL	TISSUE LEVEL	
	Ø 5.0 mm SIZE L BLOCK M1.6 black screw included OPTS160	PLATFORM R/4.8 mm SIZE L BLOCK M1.6 black screw included TS160	PLATFORM N/4.0 mm SIZE S BLOCK M1.6 black screw included TS160
1.5 mm	 OPBASE-C501-L		
2.5 mm	 OPBASE-C502-L	 TBASEC-R-L	 TBASEC-N-S
3.5 mm	 OPBASE-C503-L		

Hexagonal screwdriver.
*CEREC® is a brand of Sirona Dental Systems GmbH

PRE-MILLED RANGE

To produce customized abutment in dental laboratory, for cemented restorations

A comprehensive range for Axiom® BL and Axiom® TL

- PreFace® abutments with MEDENTiKA® holder compatible with following CAM softwares: DS CAM*, Hyperdent*, Imes-Icore*, Picasoft*, Sum 3D*, VHF® and WorkNC®.
- TI-Forms®-Abutments with Amann Girrbach® holder compatible with Straumann® CARES® M Series.
- PreFace® abutments available in 2 diameters to cover all indications.



PreFace® abutments in titanium - MEDENTiKA® holder



Removable with the gripper INEXPS/L or INEXPHS/L

Definitive titanium M1. 6 included	BONE LEVEL	TISSUE LEVEL	
		PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Ø 11.5 mm	 OPPM-M-TI-11	 TPM-M-TI-R11	 TPM-M-TI-N11
Ø 16.0 mm	 OPPM-M-TI-16	 TPM-M-TI-R16	 TPM-M-TI-N16



TI-Forms®-Abutments in titanium for M-Series - Amann Girrbach® holder for Ceramill

Definitive titanium M1. 6 included	BONE LEVEL	TISSUE LEVEL	
		PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Ø 11.5 mm	 OPPM-C-TI-11	 TPM-C-TI-R11	 TPM-C-TI-N11

Hexagonal screwdriver.
For impression components, please refer to pages 28.


*All trademarks and company names are the property of their respective owners.

MULTI-UNIT RANGE


For plural screw-retained restorations

A solution for each indication


- 2 platforms: 4.8 mm and 4.0 mm
- Secondary components common to Axiom® BL and Axiom® TL
- Easy handling

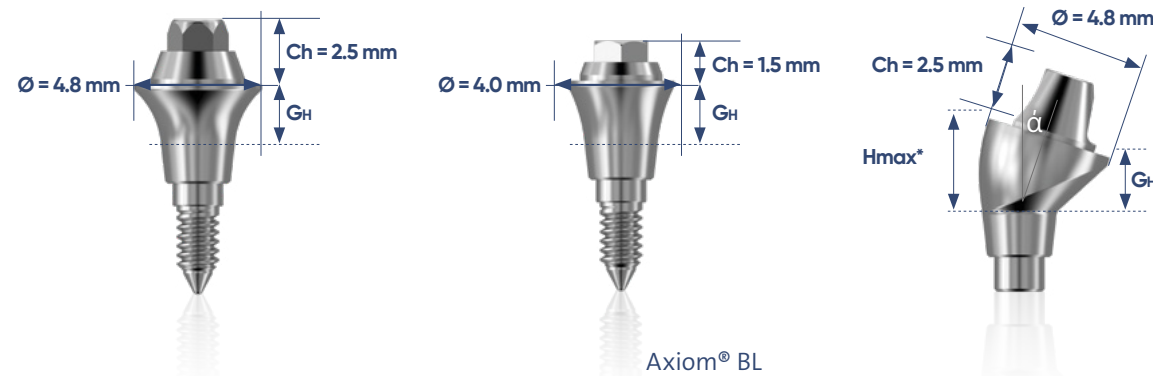

- Straight and angled abutments: 18° or 30°
- 5 gingival heights, depending on platforms
- The 20° morse taper on the platform allows greater correction of axial malalignment

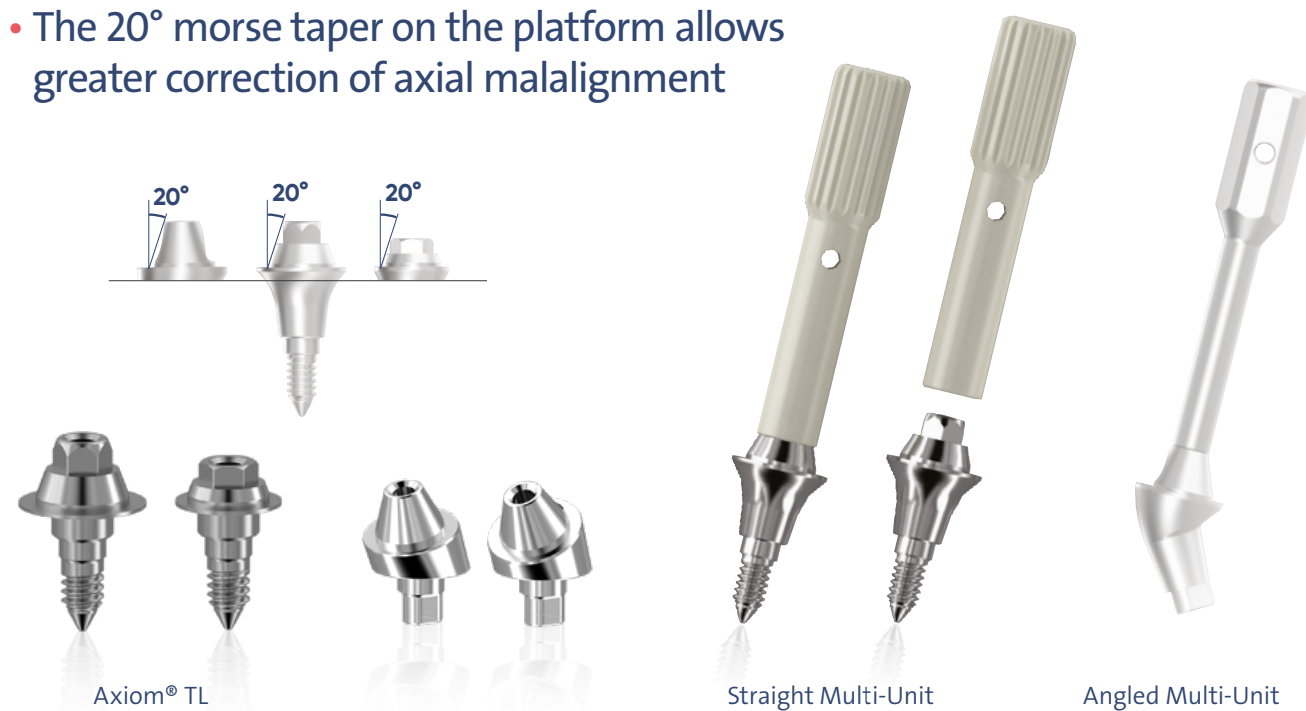

All Multi-Unit abutments are delivered with an attached holder.

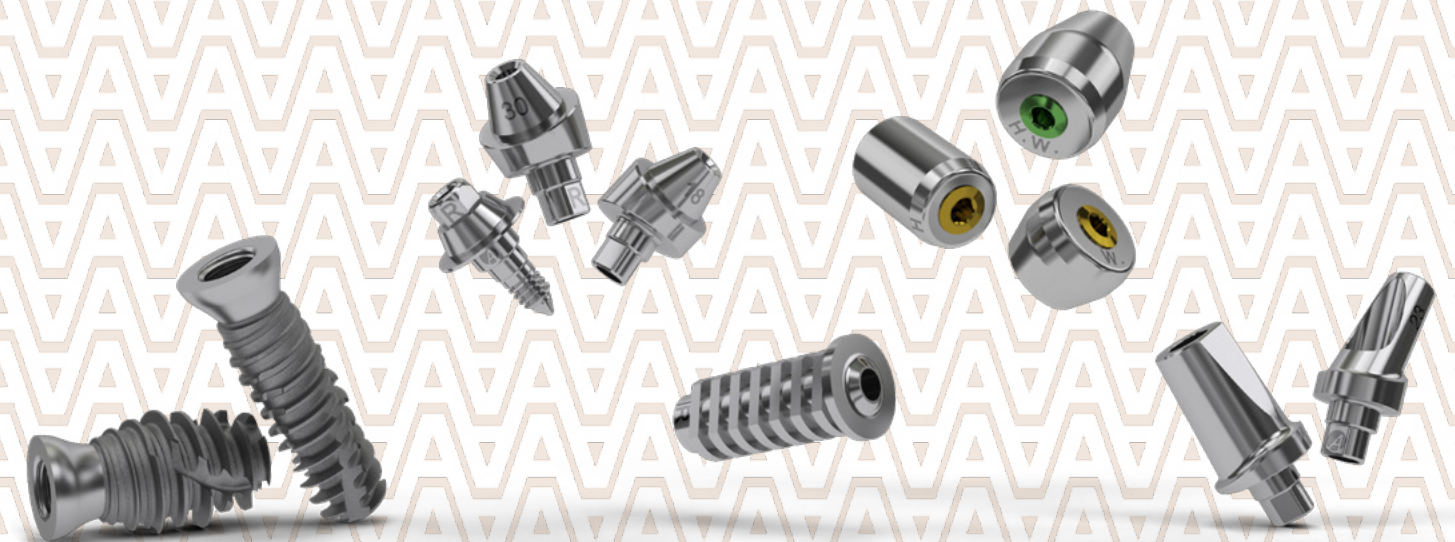

25 N.cm STERILE		BONE LEVEL				TISSUE LEVEL	
		PLATFORM R/4.8 mm		PLATFORM N/ 4.0 mm		PLATFORM R/4.8 mm	PLATFORM N/ 4.0 mm
CORONAL HEIGHT		2.5 mm		1.5 mm		2.5 mm	1.5 mm
ANGULATION	GINGIVAL HEIGHT						
STRAIGHT	0.75 mm	 OPMU0-0		 OPMUN0-0		abt Ø: 4.8  TMU0-R	abt Ø: 4.0  TMU0-N
	1.5 mm	 OPMU0-1		 OPMUN0-1			
		 OPMU0-2		 OPMUN0-2			
	2.5 mm	 OPMU0-3		 OPMUN0-3			
		 OPMU0-4		 OPMUN0-4			
	3.5 mm						
	4.5 mm						
	INDEXATION		indexed	not indexed	indexed		
18° ANGLED M1.6 black screw included Angled Multi-Unit abutment for Bone Level are removable with INEXPHS/L	1.5 mm	 OPMU18-1-IN	 OPMU18-1	-	-	abt Ø: 4.8  TMU18-R-IN	abt Ø: 4.0  TMU18-N-IN
	2.5 mm	 OPMU18-2-IN	 OPMU18-2	-	-		
	3.5 mm	 OPMU18-3-IN	 OPMU18-3	-	-		
30° ANGLED M1.6 black screw included	0.75 mm	 OPMU30-0-IN	 OPMU30-0	-	-	abt Ø: 4.8  TMU30-R-IN	abt Ø: 4.0  TMU30-N-IN
	1.5 mm	 OPMU30-1-IN	 OPMU30-1	-	-		
	2.5 mm	 OPMU30-2-IN	 OPMU30-2	-	-		
Wrench for straight Multi-Unit abutments		 MUW100		Mandrel for straight Multi-Unit abutments		 MUM100	
				Long mandrel for straight Multi-Unit abutments		 MUM100L	

 Hexagonal screwdriver.

TOP REFERENCES ARE SHOWN IN CORAL

Conventional impression on Multi-Unit

	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Multi-Unit pick-up transfer coping 1 short and 1 long screw MU included	 MUT100 MUT100-4 (X4)	 MUNT100 MUNT100-4 (X4)
Multi-Unit pop-in transfer coping Screw included	 MUT200 MUT200-4 (X4)	 MUNT200 MUNT200-4 (X4)
Multi-Unit abutment analog MU M1.4 blue screw included	 MUA100 MUA100-4 (X4)	 MUNA100 MUNA100-4(X4)

Healing on Multi-Unit

Multi-Unit protective cap (titanium) 	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Low	Ch: 4 mm  MUCAP MUCAP-4 (X4)	Ch: 3 mm  MUNCAP
High	Ch: 7 mm  MUCAP-H	Ch: 7 mm  MUNCAP-H
High & Wide	Ch: 7 mm Width: 6 mm  MUCAP-HW	Ch: 7 mm Width: 6 mm  MUNCAP-HW
Wide	Ch: 4 mm Width: 6 mm  MUCAP-W	Ch: 4 mm Width: 6 mm  MUNCAP-W

Digital impression on Multi-Unit

	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Multi-Unit digital transfer Screw included	 151-03-DT-MU	 151-04-DT-MUN
Multi-Unit metal digital transfer Screw included	 151-03-MDT	 151-04-MDT
Analog for print model Multi-Unit Laboratory screw included	 151-03-PA 151-03-PA-4 (X4)	 151-04-PA 151-04-PA-4 (X4)

Scan accessories

	MULTI-UNIT R/4.8 mm	MULTI-UNIT N/4.0 mm
Scan-Adapter Laboratory screw included  Scan-Adapter screwdriver	 151-03-SAO	 151-04-SAO
Scan-Adapter screwdriver	 SATOOL-01	

For Scan Adapter Kits, see page 106

Multi-Unit temporary components

	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Multi-Unit Titanium temporary coping MU Titanium M1.4 screw included 	 MUC100	 MUNC100
Multi-Unit PEEK temporary coping MU Titanium M1.4 screw included 	 MUC200	 MUNC200

Multi-Unit castable copings

	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Multi-Unit castable coping MU black screw MU140 included	 MUC300	 MUNC300
Multi-Unit CoCr overcast coping MU black screw MU140 included	 MUC400	 MUNC400

 Hexagonal screwdriver.

Multi-Unit temporary abutments with Angulated Access

	MULTI-UNIT Ø 4.8 mm	MULTI-UNIT Ø 4.0 mm
 M1.4 Titanium screw with AA included MUAA141 ANGULATION		
0°	 MUCAA100	 MUNCAA100
10°	 MUCAA110	 MUNCAA110
15°	 MUCAA115	 MUNCAA115
25°	 MUCAA125	 MUNCAA125

 Ball screwdriver.

Accessories associated to Angulated Access

M1.4 titanium screw with AA	 MUA141	
Blue laboratory screw M1.4 with AA	 MUA142-4 (X4)	
AA screw gripper	 AATool-M	 AATool AATool-4 (X4)
Channel rod Ø 2.5 mm for MUAA temporary or permanent prosthesis	 CR25-4 (X4)	

PACIFIC components for Multi-Unit abutments

	PLATFORM N/4.0 mm
Narrow Multi-Unit PACIFIC kit	- KITMUNPAC
Narrow Multi-Unit PACIFIC ring	 MUNPAC100
MU black screw MU140 included	
Narrow Multi-Unit PACIFIC analog	 MUNPAC110
MU laboratory screw MU141 included	
Narrow Multi-Unit PACIFIC castable coping	 MUNPAC120

Laboratory components for Multi-Unit

M 1.4 Titanium laboratory screw	 MU141	
Multi-Unit laboratory screw	 Short MUT101	 Long MUT102
	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Protective analog	 MUA200 MUA200-4 (X4)	 MUNA200
MU M1.4 blue screw (not included)		



X-Base® for plural screw-retained restorations on Multi-Unit abutments

	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
X-Base® Multi-Unit - Straight	Coronal heigh: 5 mm  MUFLEX-5		Coronal heigh: 4 mm  MUNFLEX-4	
Cast coping	 MUFLEXC		 MUNFLEXC	
Definitive screw	 MU140			
Titanium screw	 MU140Z			
Laboratory screw	Short  MU141	Long  MUT101	Longer  MUT102	
	PLATFORM R/4.8 mm		PLATFORM N/4.0 mm	
X-Base® Multi-Unit – Angulated Access	Coronal heigh: 5 mm  MUFLEX-5-AA		Coronal heigh: 4 mm  MUNFLEX-4-AA	
Cast coping AA 25°	 MUFLEXC-AA25		 MUNFLEXC-AA25	
Definitive screw	 MUA140			
Short laboratory screw	 MUA142-4 (X4)			
Angulated Access screw gripper	 AATOOL AATOOL-4 (X4)		 AATOOL-M	


ESTHETIC TITANIUM ABUTMENTS RANGE

For single-unit or multiple-unit cement-retained restorations

A wide range to cover all indications

- A choice of 4 angulations (0°, 7°, 15°, 23°)
- A design of esthetic shoulders adapted for cement-retained prosthesis

For the choice of the prosthetic component diameter, please see page 146.
Discover our Try-in abutments for Axiom® BL page 148.

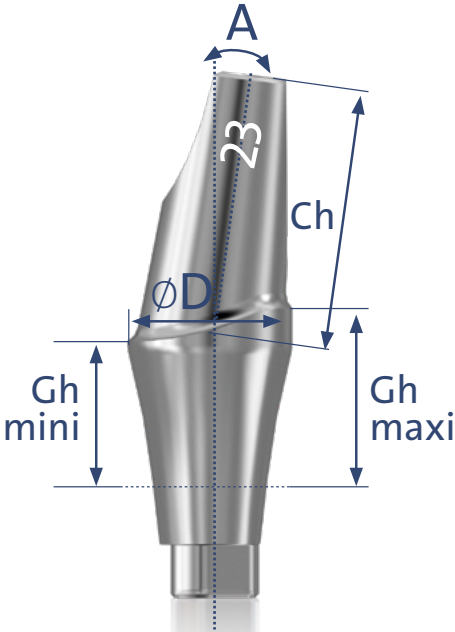

ANGULATION	GINGI- VAL HEIGHT	BONE LEVEL				TISSUE LEVEL					
		Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	Ø 4.0 mm		Ø 5.0 mm		Ø 6.0 mm	
						PLATFORM		PLATFORM		PLATFORM	
						R/4.8 mm	N/4.0 mm	R/4.8 mm	N/4.0 mm	R/4.8 mm	N/4.0 mm
M1.6 black screw included OPTS160 0°	0.75 mm	-			-						
			OPAT400	OPAT500							
	1.5 mm	-									
			OPAT410	OPAT510	OPAT610						
	2.5 mm	-									
			OPAT420	OPAT520	OPAT620						
	3.5 mm	-									
			OPAT430	OPAT530	OPAT630						
						TAT400-N		TAT500-R	TAT500-N	TAT600-R	



Removable with
the gripper
INEXPS/L
or INEXPHS/L

Hexagonal screwdriver.

ANGU- LATION	GINGI- VAL HEIGHT	BONE LEVEL				TISSUE LEVEL					
		Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	Ø 4.0 mm		Ø 5.0 mm		Ø 6.0 mm	
						PLATFORM		PLATFORM		PLATFORM	
						R/4.8mm	N/4.0mm	R/4.8mm	N/4.0mm	R/4.8mm	N/4.0mm
M1.6 black screw included OPTS160 7°	0.75 mm	-			-						
			OPAT40-7	OPAT50-7							
	1.5 mm				-						
		OPAT31-7	OPAT41-7	OPAT51-7				TAT407-N	TAT507-R	TAT507-N	
	2.5 mm				-						
		OPAT32-7	OPAT42-7	OPAT52-7							
	3.5 mm				-						
		OPAT33-7	OPAT43-7	OPAT53-7							
15°	0.75 mm	-			-						
			OPAT401	OPAT501							
	1.5 mm										
		OPAT311	OPAT411	OPAT511	OPAT611			TAT415-N	TAT515-R	TAT515-N	TAT615-R
	2.5 mm										
		OPAT321	OPAT421	OPAT521	OPAT621						
	3.5 mm										
		OPAT331	OPAT431	OPAT531	OPAT631						
23°	0.75 mm	-			-						
			OPAT402	OPAT502							
	1.5 mm	-			-						
			OPAT412	OPAT512				TAT423-N	TAT523-R	TAT523-N	
	2.5 mm	-			-						
			OPAT422	OPAT522							
	3.5 mm	-			-						
			OPAT432	OPAT532							


Note : 4.5 mm gingival height esthetic abutments are shown at page 50.

PROSTHESIS

PARTICULAR CASE:

Gingival height 4.5 mm esthetic titanium abutments have a specific gingival profile to be used in the following situations:

- Gingival height of 4.5 mm and coronal height of 7 mm
- Gingival height of 2.5 mm and coronal height of 9 mm. In order of obtain this abutment with a coronal height of 9 mm, the laboratory should rework the 2 mm cylindrical part

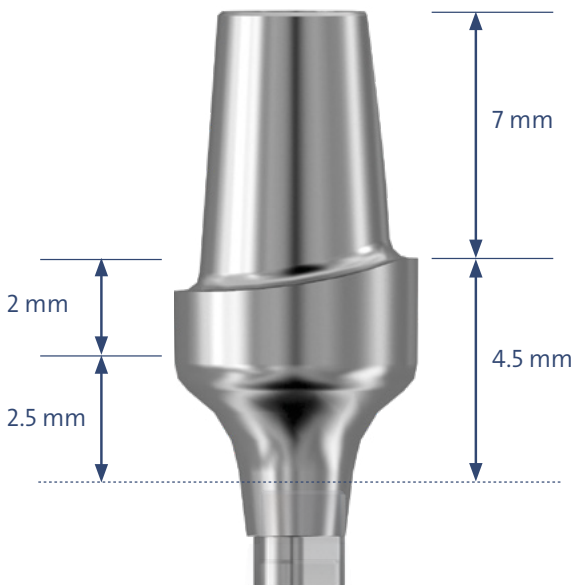

For the choice of the prosthetic component diameter, please see page 146.
Discover our Try-in abutments for Axiom® BL page 147.



M1.6 black screw included
OPTS160

		BONE LEVEL			
ANGULATION	GINGIVAL HEIGHT	Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm
0°	4.5 mm	-	 OPAT440	 OPAT540	 OPAT640
7°	4.5 mm	 OPAT34-7	 OPAT44-7	 OPAT54-7	-
15°	4.5 mm	 OPAT341	 OPAT441	 OPAT541	 OPAT641
23°	4.5 mm	-	 OPAT442	 OPAT542	-




STANDARD TITANIUM ABUTMENTS FOR AXIOM® BL

For single-unit* or multiple-unit cement-retained restorations

A wide range to cover all indications

- A choice of 3 angulations (0°, 15°, 23°)
- 3 gingival heights – 2 coronary heights
- Delivered sterile for immediate placement

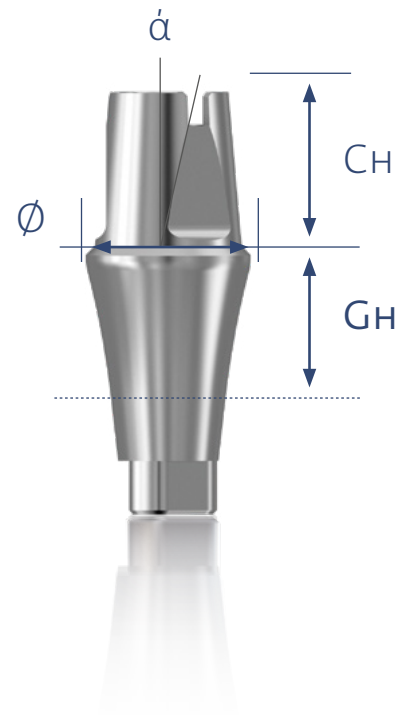

Standard Titanium abutments

ANGULATION	GINGIVAL HEIGHT	CORONAL HEIGHT	BONE LEVEL			
			Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm
M1.6 black screw included OPTS160 0° 	1.5 mm	4 mm	OPST314	OPST414	OPST514	OPST614
		6 mm	OPST316	OPST416	OPST516	OPST616
	2.5 mm	4 mm	OPST324	OPST424	OPST524	OPST624
		6 mm	OPST326	OPST426	OPST526	OPST626
	3.5 mm	4 mm	OPST334	OPST434	OPST534	OPST634
		6 mm	OPST336	OPST436	OPST536	OPST636
15° without indexation 	1.5 mm	6 mm	-	OPST416-15	OPST516-15	-
			-	OPST426-15	OPST526-15	-
	2.5 mm	6 mm	-	OPST436-15	OPST536-15	-
23° without indexation 	1.5 mm	6 mm	-	OPST416-23	OPST516-23	-
	2.5 mm	6 mm	-	OPST426-23	OPST526-23	-
	3.5 mm	6 mm	-	OPST436-23	OPST536-23	-

For the choice of the prosthetic component diameter, please see page 146.
Discover our Try-in abutments for Axiom® BL page 147.


 Hexagonal screwdriver.



Removable with the gripper
INEXPS/L or INEXPHS/L

Secondary components for Axiom® BL standard abutments

	BONE LEVEL					
	CORONAL HEIGHT	Ø 3.4 mm	Ø 4.0 mm	Ø 5.0 mm	Ø 6.0 mm	
Standard abutment transfer	 OPTT100 OPTT100-5 (X5)					
Standard abutment analog	 4 mm	OPSA304	OPSA404	OPSA504	OPSA604	
	 6 mm	OPSA306	OPSA406	OPSA506	OPSA606	
PEEK protective cap	 4 mm	OPPC304	OPPC404	OPPC504	OPPC604	
	 6 mm	OPPC306	OPPC406	OPPC506	OPPC606	
PMMA castable coping	 4 mm	Single unit anti-rotational:	OPCA304	OPCA404	OPCA504	OPCA604
		Multiple-unit rotational:	OPCR304	OPCR404	OPCR504	OPCR604
	 6 mm	Single unit anti-rotational:	OPCA306	OPCA406	OPCA506	OPCA606
		Multiple-unit rotational:	OPCR306	OPCR406	OPCR506	OPCR606

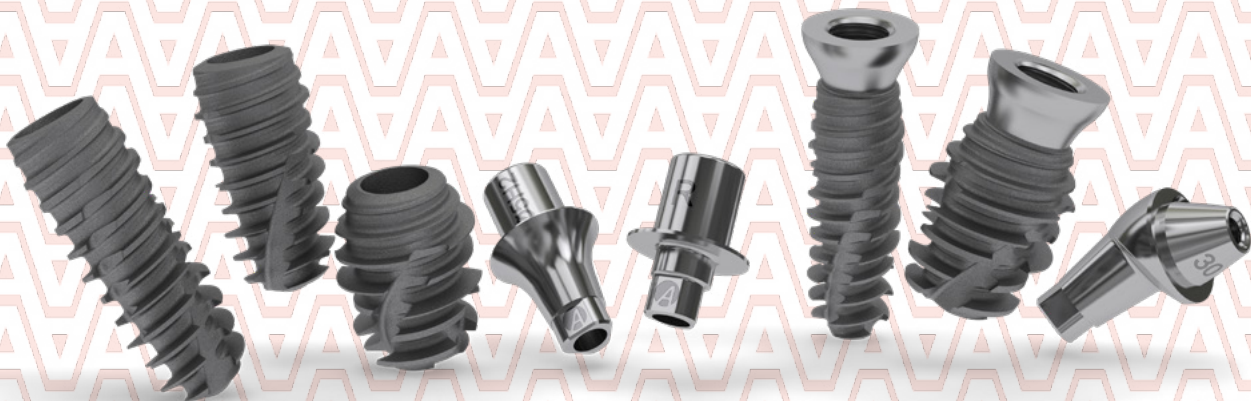

NOVALOC®


Removable overdenture retention system for edentulous patients

- 5 gingival heights, for Axiom® BL and Axiom® TL implants
- Combinaison of PEEK material and DLC* carbon-based coating offering long term performance and outstanding wear resistance
- Innovative connection system and unique convenient component handling
- Implant divergence compensation up to 40°
- Less hassle, for advanced patient comfort

* DLC = Diamond-like-Carbon



Novaloc® abutments

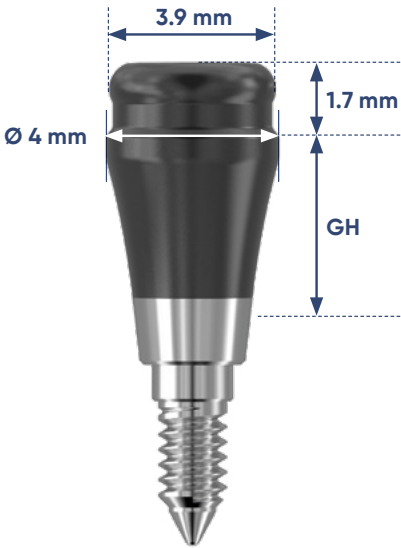

For Axiom® BL



For Axiom® TL

BONE LEVEL			TISSUE LEVEL				
	GINGIVAL HEIGHT				GINGIVAL HEIGHT	N/4.0 mm PLATFORM	R/4.8 mm PLATFORM
	1.5 mm	OPNOVA015			0.0 mm	TOPNOVA000-N	TOPNOVA000-R
	2.5 mm	OPNOVA025			1.0 mm	TOPNOVA010-N	TOPNOVA010-R
	3.5 mm	OPNOVA035			2.0 mm	TOPNOVA020-N	TOPNOVA020-R
	4.5 mm	OPNOVA045			3.0 mm	TOPNOVA030-N	TOPNOVA030-R
	5.5 mm	OPNOVA055			4.0 mm	TOPNOVA040-N	TOPNOVA040-R

NOVALOC®


Impression taking

Novaloc® impression coping	
PEEK	2010.722-NOV (X4)
Scanbody	2010.726-NOV (X4) / 2010.727-NOV (X10)

Master cast - Impression

Novaloc® model analog	
	2010.721-NOV (X4)
CAD/CAM Analog	2010.728-NOV(X4) / 2010.729-NOV (X10)
Novaloc® block out spacer	
POM	2010.723-NOV (X4)
Novaloc® processing collar	
Silicone	2010.724-NOV (X10)
Novaloc® processing insert	
POM	2010.725-NOV (X4)

Novaloc® secondary parts and special accessories

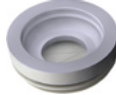Novaloc® matrix housing	Matrix: titanium	
Including mounting insert	Mounting insert: POM	2010.701-NOV (X4)
Novaloc® matrix housing - peek	Matrix: PEEK	
Including mounting insert	Mounting insert: POM	2010.702-NOV (X4)
Novaloc® matrix housing - extended	Matrix: titanium	
Including mounting insert	Mounting insert: PEEK	2010.703-NOV (X4)
Extended with attachment option		


Novaloc® inserts

Removable overdentures

Novaloc® retention insert <i>RED</i> PEEK	Extra-Light retention ≈ 300 g ≈ 0.66 lbs*	 2010.710-NOV (X4)
Novaloc® retention insert <i>WHITE</i> PEEK	Light retention ≈ 750 g ≈ 1.65 lbs*	 2010.711-NOV (X4)
Novaloc® retention insert <i>YELLOW</i> PEEK	Medium retention ≈ 1 200 g ≈ 2.65 lbs*	 2010.712-NOV (X4)
Novaloc® retention insert <i>GREEN</i> PEEK	Strong retention ≈ 1 650 g ≈ 3.64 lbs*	 2010.713-NOV (X4)
Novaloc® retention insert <i>BLUE</i> PEEK	Extra-strong retention ≈ 2 100 g ≈ 4.63 lbs*	 2010.714-NOV (X4)
Novaloc® retention insert <i>BLACK</i> PEEK	Ultra-strong retention ≈ 2 550 g ≈ 5.62 lbs*	 2010.715-NOV (X4)
Novaloc® processing package <i>TITANIUM</i> Ti / POM / PEEK / silicone	2 x 2010.701-NOV 2 x 2010.711-NOV 2 x 2010.712-NOV 2 x 2010.713-NOV 2 x 2010.724-NOV	 2010.601-NOV
Novaloc® processing package <i>PEEK</i> PEEK / POM / silicone	2 x 2010.702-NOV 2 x 2010.711-NOV 2 x 2010.712-NOV 2 x 2010.713-NOV 2 x 2010.724-NOV	 2010.611-NOV

*lbs: pounds


Novaloc® auxiliaries

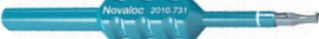Novaloc® equipment box with 3 instruments	INCLUDING 1 x 2010.731-NOV 1 x 2010.741-NOV 1 x 2010.751-NOV	 2010.101-NOV
Novaloc® processing insert removal instrument & analog repositioning aid <i>BLUE</i>		 2010.731-NOV
Novaloc® retention insert instrument (mounting and demounting) <i>BROWN</i>		 2010.741-NOV
Novaloc® matrix housing extraction instrument <i>GREY</i>		 2010.751-NOV

LOCATOR®

Removable overdenture retention system for edentulous patients

- 5 gingival heights, for Axiom® BL and Axiom® TL implants
- Proven dual retentive system offering stability, hardness and wear resistance
- Pivoting technology and self-aligning design for advanced patient comfort and usability
- Divergence compensation up to 40°
- Compatible with original Locator® and Locator® fixed attachment systems

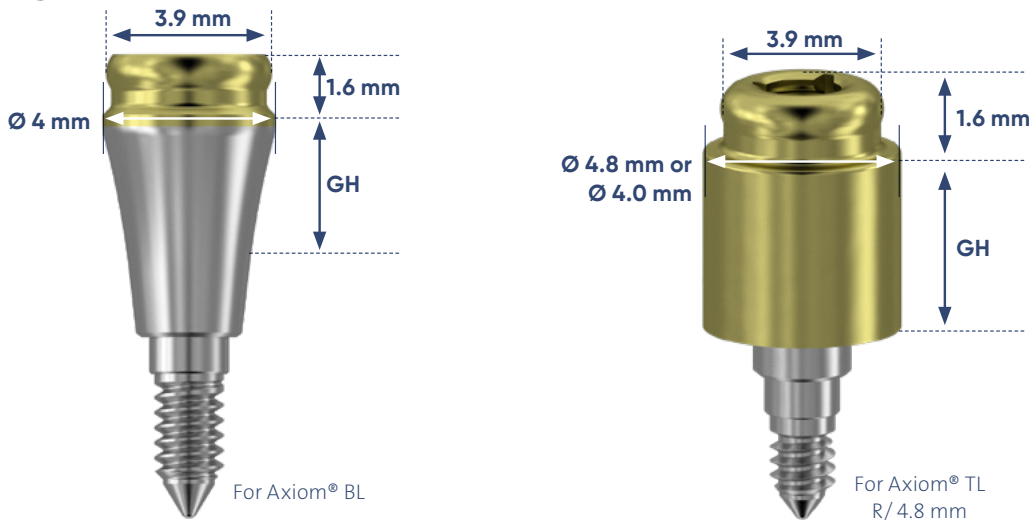

Locator® abutments

BONE LEVEL			TISSUE LEVEL			
	GINGIVAL HEIGHT			GINGIVAL HEIGHT	N/Ø 4.0 mm PLATFORM	R/Ø 4.8 mm PLATFORM
	1.5 mm	O1697				
	2.5 mm	O1698				
	3.5 mm	O1699				
	4.5 mm	O1016				
	5.5 mm	O1017				

LOCATOR®

Impression taking

Locator® impression coping	 O8505 (X4)
LOCATOR® scanbody	 O62.5057V2 (X2) O62.5057V4 (X4) O62.5057V10 (X10)

Master cast - Impression

Locator® analog	 O8530 (X4)
-----------------	----------------

LOCATOR® ATTACHMENT SYSTEM

Locator® secondary parts and special accessories

STANDARD RANGE	INSERT RETENTION	UP TO 20° DIVERGENCE COMPENSATION BETWEEN 2 IMPLANTS
Locator® male processing package Titanium, polyethylene, nylon	Standard range retention including: 1 x O8524 1 x O8527 1 x O8529 Spacer ring Black housing	 O8519-2 (X2)
Locator® replacement male insert CLEAR Nylon	Standard retention ≈ 2 268 g ≈ 5.0 lbs*	 O8524 (X4)
Locator® replacement male insert PINK Nylon	Light retention ≈ 1 361 g ≈ 3.0 lbs*	 O8527 (X4)
Locator® replacement male insert BLUE Nylon	Extra-Light retention ≈ 680 g ≈ 1.5 lbs*	 O8529 (X4)

Locator® inserts secondary parts

EXTENDED RANGE	INSERT RETENTION	UP TO 40° DIVERGENCE COMPENSATION BETWEEN 2 IMPLANTS
Locator® male processing package Titanium, polyethylene, nylon	Extended range retention including: 1 x O8558 1 x O8548 1 x O8915 Spacer ring Black housing	O8540-2 (X2)
Locator® replacement male insert GREY Nylon	Zero retention 0 g 0 lbs*	O8558 (X4)
Locator® replacement male insert RED Nylon	Extra-light retention ≈ 226 – 680 g ≈ 0.5 – 1.5 lbs*	O8548 (X4)
Locator® replacement male insert ORANGE Nylon	Light retention ≈ 907 g ≈ 2.0 lbs*	O8915 (X4)
Locator® replacement male insert GREEN Nylon	Medium retention ≈ 1361 – 1814 g ≈ 3.0 – 4.0 lbs*	O8547 (X4)

*lbs: pounds

Locator® auxiliaries

Locator® torque wrench insert driver Length: 15 mm	Stainless steel	O8260
Locator® core tool three-part Length: 100 mm	Stainless steel	O8393
Locator® torque wrench insert latch Length: 23 mm	Stainless steel	O8913

LOCATOR® FIXED ATTACHMENT SYSTEM

A financially accessible, FIXED full arch treatment solution using the same abutments and workflow as the original LOCATOR® Removable Attachment System




LOCATOR FIXED®		
LOCATOR FIXED® AUXILIARIES		NEW
LOCATOR FIXED® Insert, Four Unit (GREEN)	048.964 (x2) 048.965 (x10)	
LOCATOR FIXED® Insert, Anterior/Posterior	048.966 (x2) 048.967 (x10)	
LOCATOR FIXED® Insert, Mid Arch	048.968 (x2) 048.969 (x10)	
LOCATOR FIXED® Denture Attachment Housing	048.970 (x4) 048.971 (x10)	
LOCATOR FIXED® Processing Package, Four Unit	048.972	
LOCATOR FIXED® Processing Package, Anterior/Posterior	048.973	
LOCATOR FIXED® Processing Package, Mid Arch	048.974	

LOCATOR FIXED®		
LOCATOR FIXED® TOOLS		
LOCATOR FIXED® Seating and Removal Tool	046.828	
LOCATOR® Systems Core Tool	046.829	
LOCATOR FIXED® Seating Tip Assembly	10157-01	
LOCATOR FIXED® Removal Tip Assembly	10156-01	
LOCATOR FIXED® Removal Tip Loop	10159-01	
LOCATOR FIXED® Seating Tip Cushion	10158-01	
LOCATOR FIXED® 2.5mm Hex Wrench	10154-01	
LOCATOR FIXED® Tip Wrench	10155-01	
LOCATOR FIXED® Optional Removal Tip	10095	
LOCATOR FIXED® Storage Container	046.830	

AxIN®

Simeda® customized screw-retained tooth WITH ANGULATED ACCESS

- For all sectors*
- No glue, no sealing cement: biological safety
- Narrow angulated channel: 2 mm
- Screw integrated in the restoration: easy handling
- AxIN® base and crown delivered in one piece

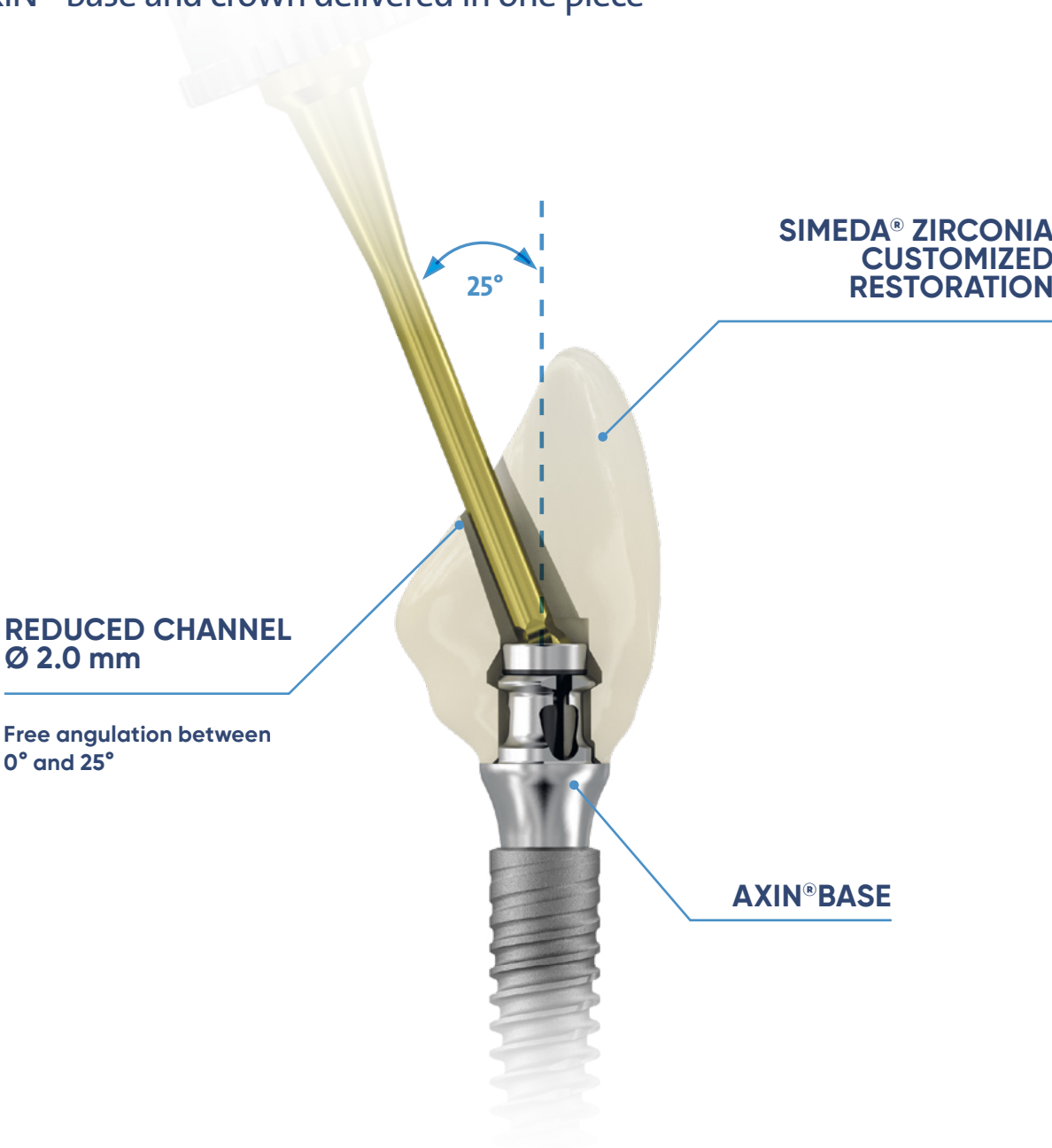

Temporary components for AxIN®


 AxIN® base	GINGIVAL HEIGHT	BONE LEVEL		TISSUE LEVEL	
		Ø 4.0 mm	Ø 5.0 mm	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
 Removable with the gripper INEXPS/L	1.5 mm Dedicated black screw included AXIN152-27-S1	 AXIN152-27B41	 AXIN152-27B51	 AXIN156-02-B	 AXIN156-01-B
	2.5 mm Dedicated black screw included AXIN152-27-S2	 AXIN152-27B42	 AXIN152-27B52	Dedicated screw included  AXIN156-0X-S	
	ANGULATION	BONE LEVEL		TISSUE LEVEL	
		Ø 4.0 mm	Ø 5.0 mm	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
AxIN® temporary crown	10° 	AXIN-C-40-10	AXIN-C-50-10	AXIN-C-48-10	AXIN-C-40-10
	15° 	AXIN-C-40-15	AXIN-C-50-15	AXIN-C-48-15	AXIN-C-40-15
	25° 	AXIN-C-40-25	AXIN-C-50-25	AXIN-C-48-25	AXIN-C-40-25
AxIN® protection cap		AXIN-PCC-40	AXIN-PCC-50	AXIN-PCC-48	AXIN-PCC-40

● Ball screwdriver.
For impression components, please refer to pages 28.

AxIN® laboratory parts

	BONE LEVEL		TISSUE LEVEL	
GINGIVAL HEIGHT	Ø 4.0 mm	Ø 5.0 mm	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
1.5 mm With laboratory screw AXIN152-27SL1 inserted	AXIN152-27L41	AXIN152-27L51	AXIN156-02-L	AXIN156-01-L
2.5 mm With laboratory screw AXIN152-27SL2 inserted	AXIN152-27L42	AXIN152-27L52	With laboratory screw AXIN156-0X-SL inserted	
Channel rod Ø 2.0 mm for temporary or perma- nent AxIN® prosthesis	CR20-4 (X4)			


● Ball screwdriver.

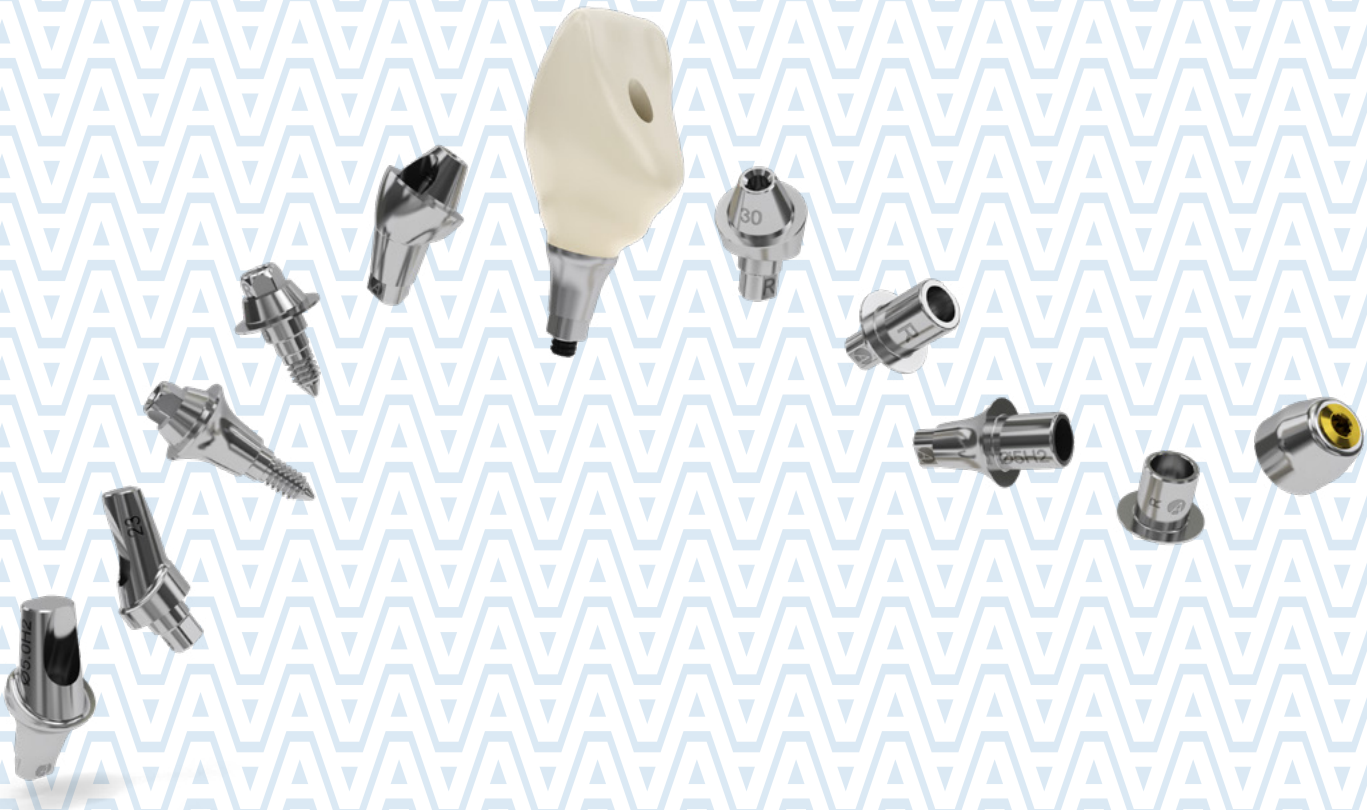

SIMEDA® CAD/CAM PLURAL PROSTHESIS

Customized screw-retained bars & bridges on Multi-Unit abutments or direct on Axiom® Tissue Level implants, straight or Angulated Access up to 25°

- Perfect fit of customized prosthesis
- Angulated Access: enhanced biosafety and screw channels outside visible or fragile areas
- No intermediate abutments on Axiom® Tissue Level: cost efficiency
- No Ti-base: time saving and esthetics
- Ease of use: delivered with an Angulated Access screw gripper
- For implant bridges and suprastructures in Titanium, Cobalt - Chromium or Zirconia

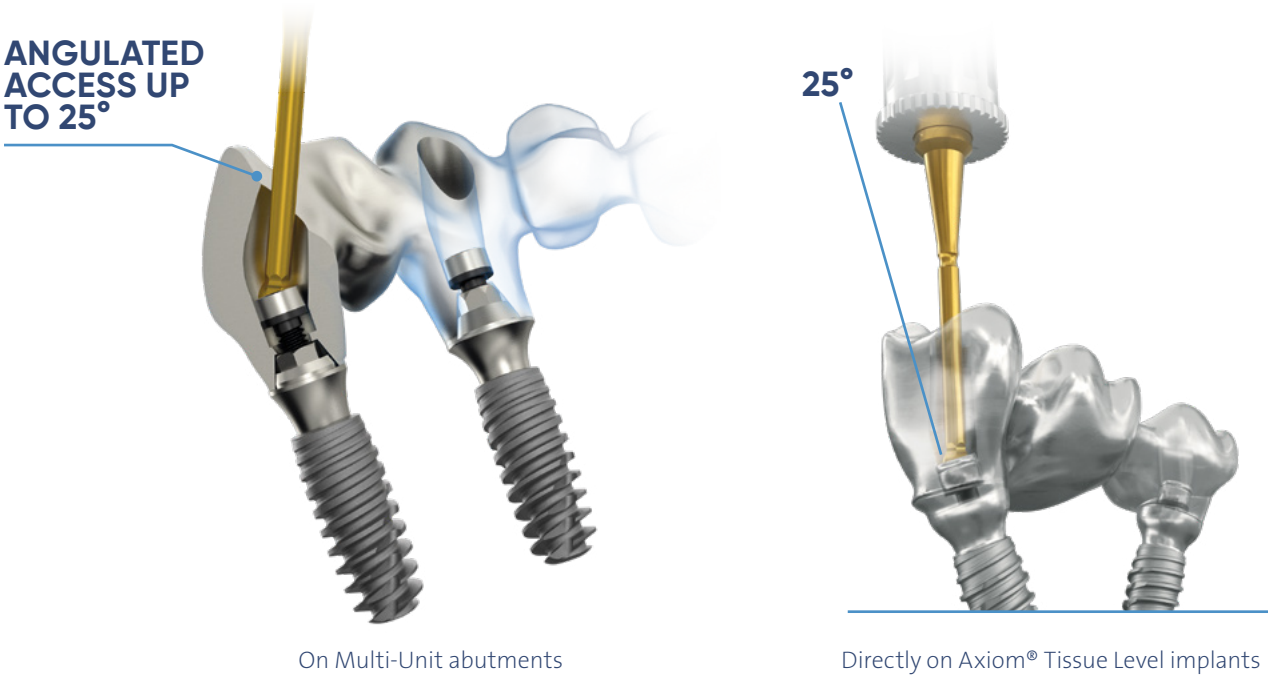

Screws for Simeda® prosthesis on Multi-Unit abutments

M1.4 titanium screw with AA	 MUAA141
Blue laboratory screw M1.4 with AA	 MUAA142-4 (X4)

Accessories for Simeda® prosthesis directly on Axiom® Tissue Level

	WITH AA	STRAIGHT ACCESS	
		SINGLE	Plural
Definitive black screw	 TSAA160P	 TS160	 TS160P
Definitive titanium screw (for zirconia prosthesis)	 TSAA161P		 TS161P
Laboratory screw	 TSAA163P-2	 TS163  TS162	 TS163P-2 (X2)  TS162P-2 (X2)

Accessories for Angulated Access prosthesis

AA screw gripper	 AATOOL AATOOL-4 (X4)	 AATOOL-M
Channel rod Ø 2.5 mm for MUAA temporary or permanent prosthesis	 CR25-4 (X4)	

SURGERY KIT

1 single compact surgical kit for all Axiom® BL and Axiom® TL, X3, X3S, REG and PX profiles


- Simplicity
- Freedom and flexibility
 - Possibility to adapt the implant choice during surgery
- User friendly
 - Complementary solutions



- Hard biocompatible coated drills
- Surface treatment optimizing depth marks visibility in mouth


Axiom® surgery kit - Full version

- Axiom Bone Level and Tissue Level
- All profiles - X3 - X3S - REG - PX
- All diameters - from Ø 3.4 to 6.4 mm

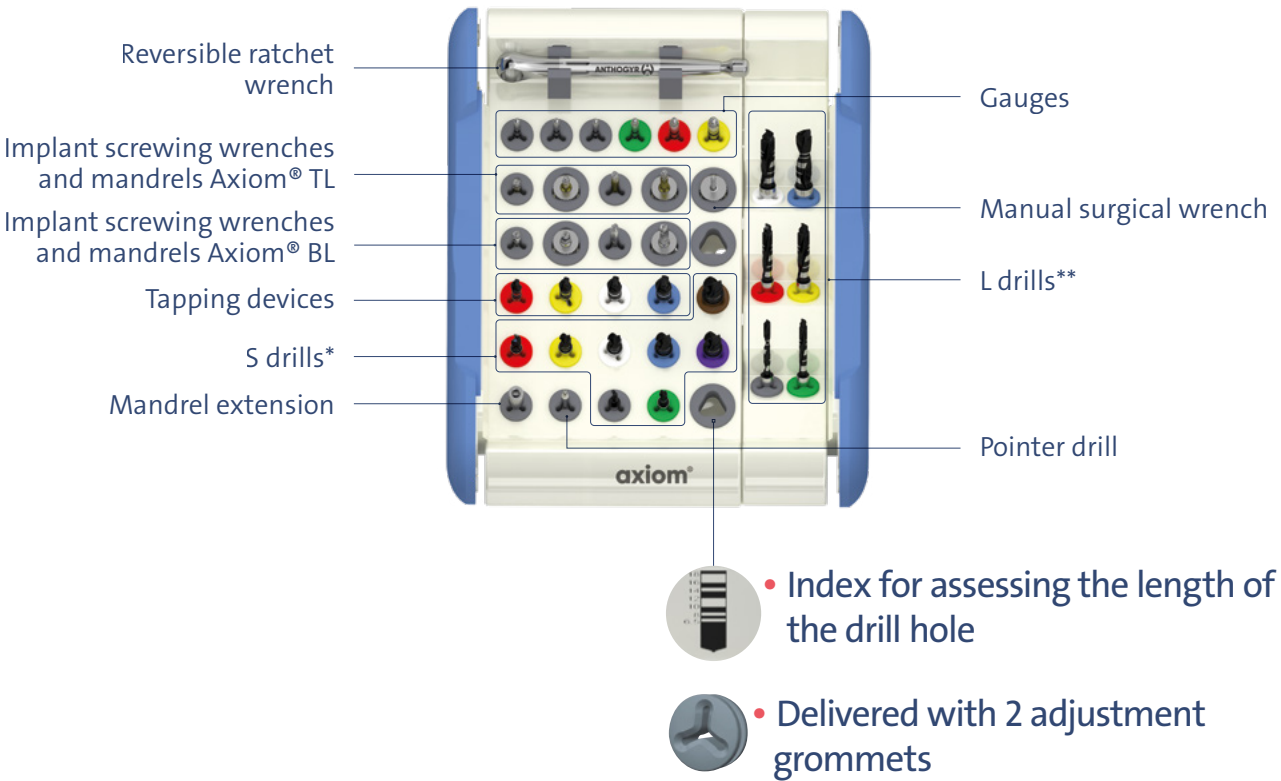

Axiom® surgery kit - Bone Level Version

- Axiom Bone Level only
- All profiles - X3 - X3S - REG - PX
- All diameters - from Ø 3.4 to 6.4 mm

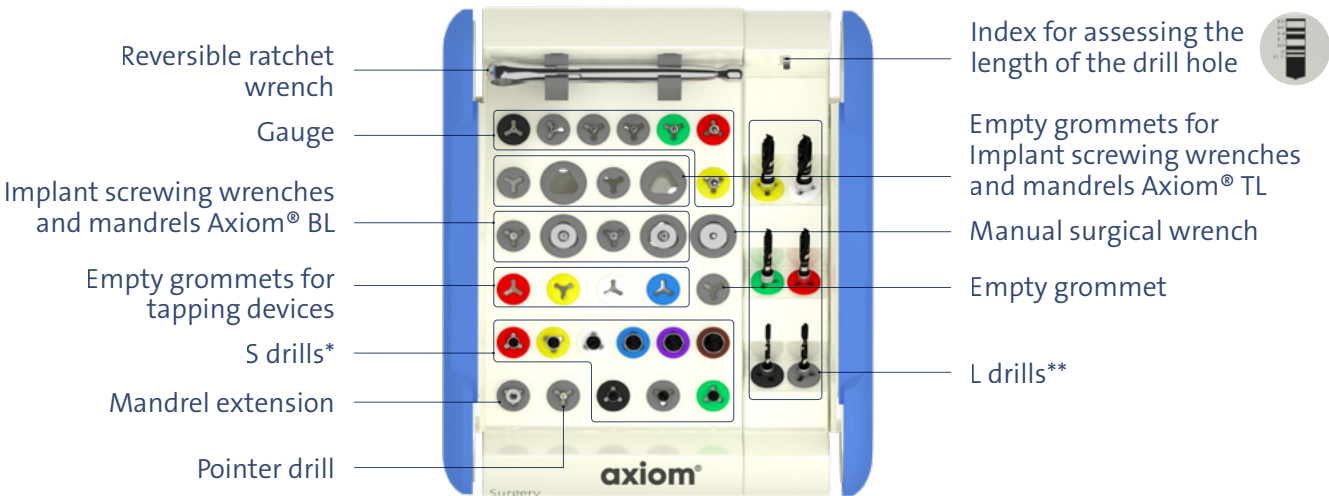

Surgery kits, instruments and accessories

Surgical kits references depending upon local availability.

		USE
NEW Axiom® surgery kit - Full version	 INMODOPS3F	For Axiom® BL and TL X3, X3S, REG and PX profiles All diameters, from Ø 3.4 to 6.4 mm
Empty Axiom® surgery kit	NEW INMODOPS3FV	
NEW Axiom® surgery kit - Bone Level version	 INMODOPS3S	For Axiom® BL only X3, X3S, REG and PX profiles All diameters, from Ø 3.4 to 6.4 mm
Axiom® surgery kit	 INMODOPS3L	For Axiom® BL and TL, X3, REG and PX profiles All diameters, from Ø 3.4 to 6.4 mm
Empty Axiom® surgery kit	INMODOPS3LV	
Manual surgical wrench	 OPCS100	
Mandrel extension	 INEXM	
Mandrel wrench	 INCPM*	
Reversible ratchet wrench	 INCC	
Surgical dynamometrical wrench*	 INCCDC*	0 to 80 N.cm
Axiom® dual function depth gauge*	 OPJC001*	
Drill guide*	 INGPPA*	
Angled drilling guide*	 INGFA*	

Universal surgical instrument	1 universal surgical instrument + 1 trilobe mandrel extension + 1 trilobe mandrel for Axiom® Bone Level implants (2.8 excluded)	 INUSI
Trilobe mandrel extension		INEXMOST
Trilobe mandrel for Axiom® Bone Level implants (2.8 excluded)		OPMVTOST
Trilobe mandrel for Axiom® Tissue Level implants		TIMTOST

	DIAMETER	SHORT	LONG
Round bur*		 INFB20*	
Pointer drill	1.5 mm	 OPPO15150	
Lindemann bur*	2.0 mm	 OPR20*	
Initial drill	1.5 mm	 OPFI15S	 OPFI15L
	2.0 mm	 OPFI20S	 OPFI20L
Step drill	2.0/ 2.4 mm ●	 OPFE24S	 OPFE24L
	2.4/ 3.0 mm ●	 OPFE30S	 OPFE30L
	3.0/ 3.6 mm ●	 OPFE36S	 OPFE36L
	3.6/ 4.2 mm ○	 OPFE42S	 OPFE42L
	4.2/ 4.8 mm ●	 OPFE48S	 OPFE48L
	4.8 / 5.4 mm ●	 OPFE54S**	-
	5.4 / 6.0 mm ●	 OPFE60S**	-
Pack of 6 step drills	Assorted	 OPFES-6	 OPFEL-6
Pack of 12 step drills (6 S + 6 L)	Assorted	 OPFESL-12	

* Not included in INMODOPS3 and INMODOPS3L kits.
** Only in the INMODOPS3L kit.

For Axiom® PX taps, see page 85

	DIAMETER	
Implant tap, for Axiom® REG	3.4 mm ●	 OPTA34L
	4.0 mm ●	 OPTA40L
	4.6 mm ○	 OPTA46L
	5.2 mm ●	 OPTA52L
Pack of 4 taps, for Axiom® REG	Assorted	 OPTA-4
Gauge	2.0 mm	 OPJD020
	2.4 mm	 OPJD024
	3.0 mm	 OPJD030
	3.6 mm	 OPJD036
	4.2 mm*	 OPJD042*
	4.8 mm*	 OPJD048*
	5.4 mm*	 OPJD054*
	6.0 mm*	 OPJD060*

* Not included in INMODOPS3 and INMODOPS3L kits.

	SHORT	MEDIUM	LONG
Axiom® BL implant mandrel	 OPMV180 23 mm	 OPMV215* 27 mm	 OPMV250 32 mm
Axiom® BL manual wrench	 OPCV060 21 mm	 OPCV110* 27 mm	 OPCV160 31 mm
Axiom® BL implant manual screw wrench (with screwdriver)		 INPIM*	 INPIL*
Axiom® TL implant mandrel	 TIM100S 19 mm		 TIM100L 29 mm
Axiom® TL abutment wrench	 TIW100S 22 mm		 TIW100L 32 mm
Screw extraction mandrel*	 INMDS*		 INMDL*
Pack of Axiom® TL mandrels and wrenches	 TIMW100-4		

* Not included in INMODOPS3 and INMODOPS3L kits.

AXIOM® IMPLANT PRE-POSITIONING MG KIT

To EASILY mark implant position regarding an adjacent tooth or an implant

- Quick and easy guidance
- Immediate solution
- Reusable and cost efficient
- Time saving

CENTERING RINGS

- To position referring to an adjacent tooth
- 2 diameters


POINTER DRILL

PARALLELING DRILL GUIDE

SMART REVERSIBLE 2-IN-1 ROD

- To position referring to an implant bed
- To position referring to an Axiom® BL implant

Axiom® MG Kit and instruments

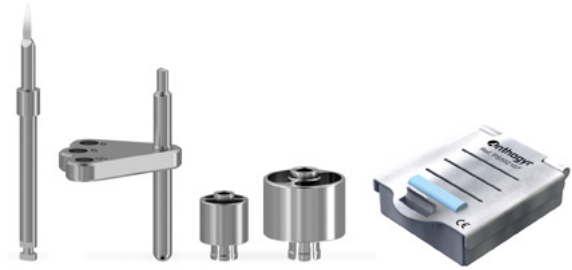Axiom® implant pre-positioning MG Kit	 INGUIDE
Parallelizing drill guide with smart 2-in-1 rod	 INGPPA
Pointer drill Ø1.5	 OPPO15
Centering ring Ø8.0	 OPR8
Centering ring Ø10	 OPR10
Empty micro-cassette for surgical instrument storage and sterilization	 PS502107

DENSE BONE KIT

To facilitate the placement of Axiom® implants in high bone densities

- Compact kit dedicated to dense bone indications, to be used in addition to Axiom® surgery kit
- Short step cortical drills to reduce cortical bone constraint, for Axiom® implants
- Specific Axiom® PX taps for D1 bone



Axiom® dense Bone kit and instruments


Axiom® Dense Bone surgery kit	 INMODCP
Empty Axiom® Dense Bone surgery kit	 INMODCPV

Axiom® cortical step drill Ø 3.0 / 3.2 S ●	 PXFE32S
Axiom® cortical step drill Ø 3.6 / 3.8 S ●	 PXFE38S
Axiom® cortical step drill Ø 4.2 / 4.4 S ⊖	 PXFE44S
Axiom® cortical step drill Ø 4.8 / 5.0 S ⊖	 PXFE50S
Pack of 4 Axiom® cortical step drills S	 PXFES-4
Axiom® PX tap Ø 3.4 L for dense bone ●	 PXTA34L
Axiom® PX tap Ø 4.0 L for dense bone ●	 PXTA40L
Axiom® PX tap Ø 4.6 L for dense bone ⊖	 PXTA46L
Axiom® PX tap Ø 5.2 L for dense bone ⊖	 PXTA52L
Pack of 4 Axiom® PX taps L for dense bone	 PXTAL-4


COUNTERSINKS

To optimize bone crest preparation for easier implant and prosthetic placement

- Invasiveness management
- Accurate and safe guidance
- Guarantee of a proper final seating

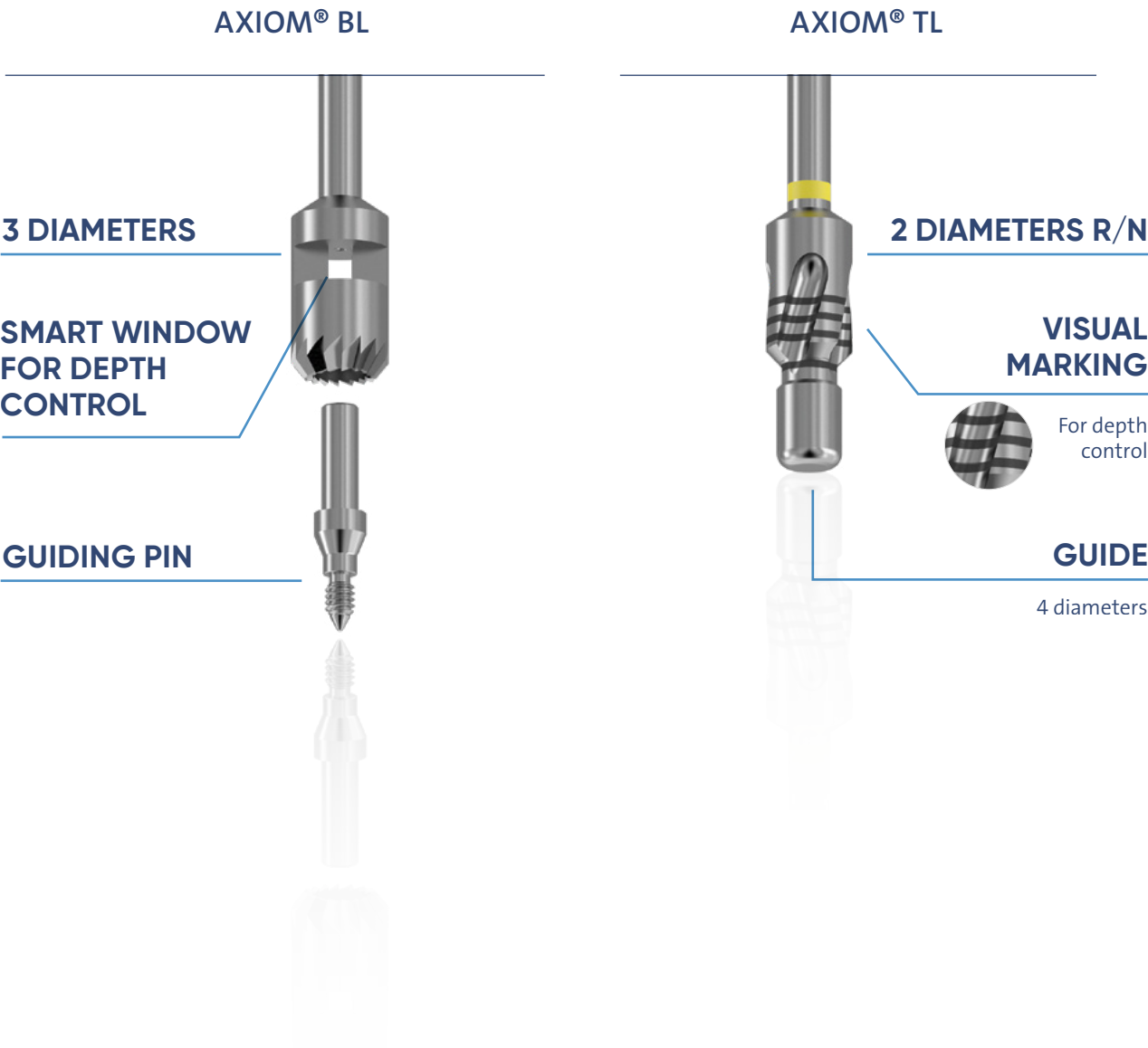

Axiom® countersinks

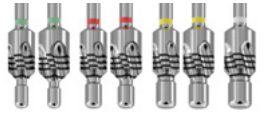	DIAMETER	BONE LEVEL	TISSUE LEVEL	
			PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Pack of 3 Axiom® BL countersinks with pins	Assorted	 OPFF-3	-	-
Axiom® BL countersink with pin	4.5 mm	 OPFF45	-	-
	5.3 mm	 OPFF53	-	-
	6.6 mm	 OPFF66	-	-
Axiom® BL countersink pin 	-	 OPFFP	-	-
Pack of 7 Axiom® TL countersinks	Assorted	-	 TOPFF-7	
Axiom® TL countersink	2.4 mm ●	-	 TOPFF24R	 TOPFF24N
	3.0 mm ●	-	 TOPFF30R	 TOPFF30N
	3.6 mm ●	-	 TOPFF36R	 TOPFF36N
	4.2 mm ○	-	 TOPFF42R	-

DRILL STOP KIT

Ergonomic and user-friendly

- Easy drill stop pick-up with contra-angle
- Anti-falling system
- Color coding facilitating identification
- For all Axiom® Bone Level and Axiom® Tissue Level implants




DRILL STOP KIT	AXIOM® MULTI LEVEL®*
Axiom® drill stop kit For all Axiom® diameters, from 3.4 to 6.4 mm 18 short drill stops (S) and 24 long drill stops (L)	 INKITOPDSL
Empty Axiom® drill stop box For 42 drill stops	INKITOPDSL
Axiom® drill stop kit For diameters 3.4 to 5.2 mm 12 short drill stops (S) and 24 long drill stops (L)	 INKITOPDS
Empty Axiom® drill stop box for 36 drill stops	INKITOPDS

	STOP FOR S DRILLS						STOP FOR L DRILLS			
Drill diameter										
Implant length	Ø 2.0 / 2.4 / 3.0 mm	Ø 3.6 mm	Ø 4.2 mm	Ø 4.8 mm	Ø 5.4 mm	Ø 6.0 mm	Ø 2.0 / 2.4 / 3.0 mm	Ø 3.6 mm	Ø 4.2 mm	Ø 4.8 mm
6.5 mm	OPB3006C Ø 3.0 mm	OPB3606C Ø 3.6 mm	OPB4206C Ø 4.2 mm	OPB4806C Ø 4.8 mm	OPB5406C Ø 5.4 mm	OPB6006C Ø 6.0 mm	OPB3006L Ø 3.0 mm	OPB3606L Ø 3.6 mm	OPB4206L Ø 4.2 mm	OPB4806L Ø 4.8 mm
8 mm	OPB3008C Ø 3.0 mm	OPB3608C Ø 3.6 mm	OPB4208C Ø 4.2 mm	OPB4808C Ø 4.8 mm	OPB5408C Ø 5.4 mm	OPB6008C Ø 6.0 mm	OPB3008L Ø 3.0 mm	OPB3608L Ø 3.6 mm	OPB4208L Ø 4.2 mm	OPB4808L Ø 4.8 mm
10 mm	OPB3010C Ø 3.0 mm	OPB3610C Ø 3.6 mm	OPB4210C Ø 4.2 mm	OPB4810C Ø 4.8 mm	OPB5410C Ø 5.4 mm	OPB6010C Ø 6.0 mm	OPB3010L Ø 3.0 mm	OPB3610L Ø 3.6 mm	OPB4210L Ø 4.2 mm	OPB4810L Ø 4.8 mm
12 mm							OPB3012L Ø 3.0 mm	OPB3612L Ø 3.6 mm	OPB4212L Ø 4.2 mm	OPB4812L Ø 4.8 mm
14 mm							OPB3014L Ø 3.0 mm	OPB3614L Ø 3.6 mm	OPB4214L Ø 4.2 mm	OPB4814L Ø 4.8 mm
16 mm							OPB3016L Ø 3.0 mm	OPB3616L Ø 3.6 mm		
18 mm							OPB3018L Ø 3.0 mm	OPB3618L Ø 3.6 mm		

INTEGRAL GUIDED SURGERY


Fully guided surgery solution to cover main indications with Axiom® Bone and Tissue Level implants, X3, REG and PX profiles

- Ergonomic and efficient system
- Accurate components adjustment for high precision guidance
- Follow the path with intuitive color coding
- Open compatibility with major planning software

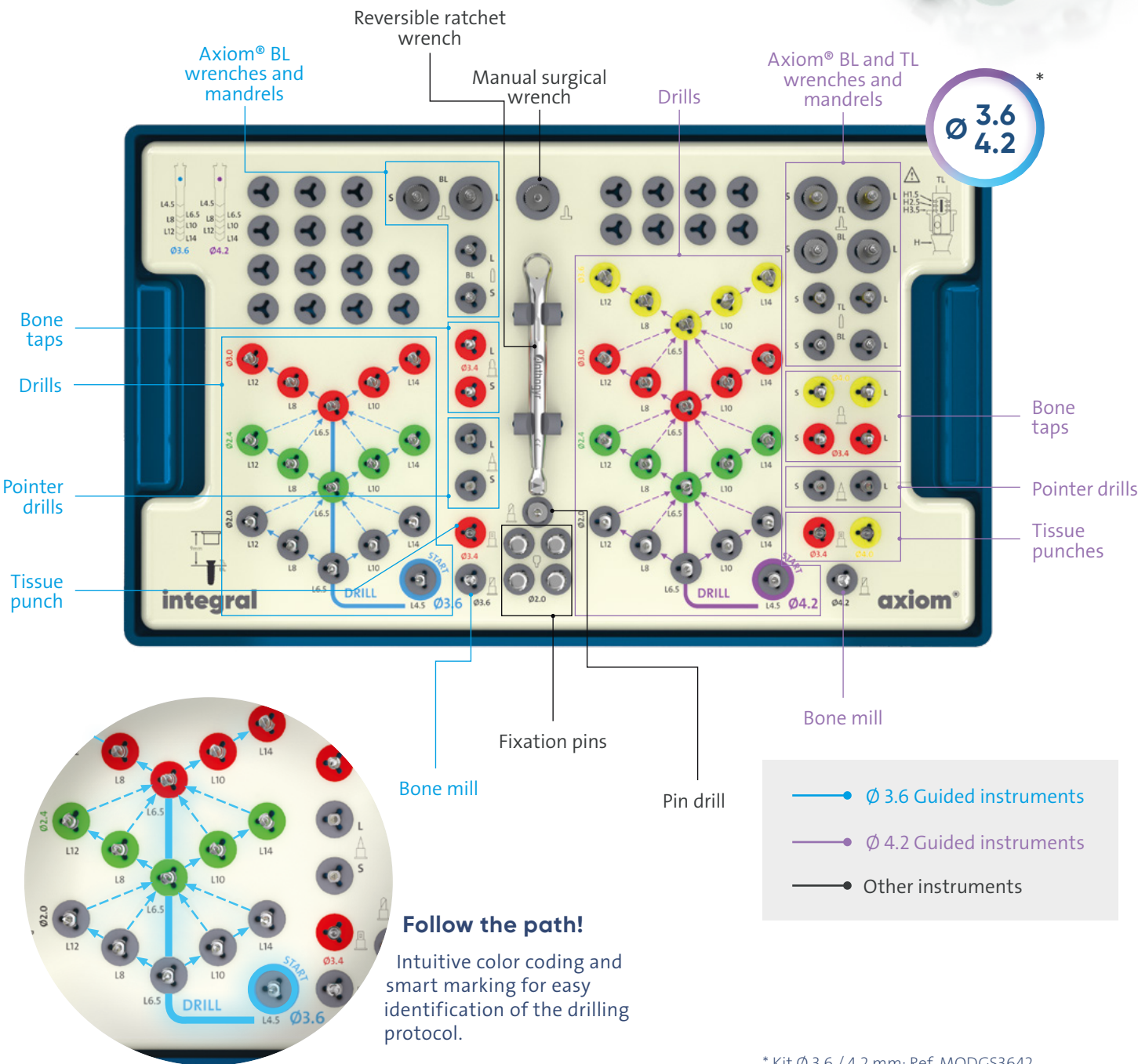


INTEGRAL SURGICAL KITS

INTEGRAL Guided surgery kit, Ø 3.6 and 4.2

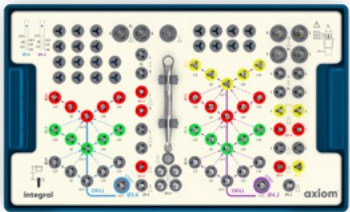


MODGS3642

INTEGRAL Guided surgery kit, Ø 3.6

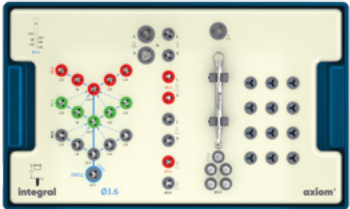


MODGS36

INTEGRAL Guided surgery kit, Ø 4.2

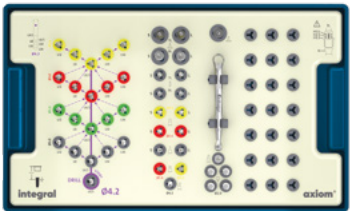

MODGS42

INTEGRAL Guided surgery kit, Ø 5.0

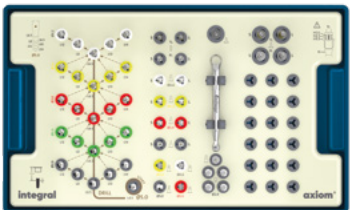

MODGS50

Empty INTEGRAL Guided surgery kit Ø 3.6 and Ø 4.2

MODGS3642V

Empty INTEGRAL Guided surgery kit Ø 3.6

MODGS36V

Empty INTEGRAL Guided surgery kit Ø 4.2

MODGS42V

Empty INTEGRAL Guided surgery kit Ø 5.0

MODGS50V

Guided surgery sleeve insertion tool




OPGSTOOL-01

INTEGRAL GUIDE SLEEVES				
		Ø 3.6 mm ●	Ø 4.2 mm ●	Ø 5.0 mm ●
Guide sleeves *		 OPGS36 OPGS36-4 (X4) OPGS36-10 (X10)	 OPGS42 OPGS42-4 (X4) OPGS42-10 (X10)	 OPGS50 OPGS50-4 (X4) OPGS50-10 (X10)
INTEGRAL INSTRUMENTS				
GUIDED TISSUE PUNCHES		For sleeves Ø 3.6 mm ●	For sleeves Ø 4.2 mm ●	For sleeves Ø 5.0 mm ●
Guided tissue punches	Ø 3.4 mm	 GTP34G36S	 GTP34G42S	 GTP34G50S
	Ø 4.0 mm	-	 GTP40G42S	 GTP40G50S
	Ø 4.6 mm	-	-	 GTP46G50S
GUIDED BURS		For sleeves Ø 3.6 mm ●	For sleeves Ø 4.2 mm ●	For sleeves Ø 5.0 mm ●
Guided bone mills		 GBAG36S	 GBAG42S	 GBAG50S
Guided cortical burs Ø 2.4 mm*	Short	 GCBG36S	 GCBG42S	 GCBG50S
	Long	 GCBG36L	 GCBG42L	 GCBG50L
GUIDED POINTER DRILLS		For sleeves Ø 3.6 mm ●	For sleeves Ø 4.2 mm ●	For sleeves Ø 5.0 mm ●
Guided pointer drills	Short	 GPDG36S	 GPDG42S	 GPDG50S
	Long	 GPDG36L	 GPDG42L	 GPDG50L
GUIDED INITIAL DRILLS		For sleeves Ø 3.6 mm ●	For sleeves Ø 4.2 mm ●	For sleeves Ø 5.0 mm ●
Guided initial drills Ø 2.0 mm	Length 4.5 mm	 GFI20L4G36	 GFI20L4G42	 GFI20L4G50
	Length 6.5 mm	 GFI20L6G36	 GFI20L6G42	 GFI20L6G50
	Length 8 mm	 GFI20L8G36	 GFI20L8G42	 GFI20L8G50
	Length 10 mm	 GFI20L10G36	 GFI20L10G42	 GFI20L10G50
	Length 12 mm	 GFI20L12G36	 GFI20L12G42	 GFI20L12G50
	Length 14 mm	 GFI20L14G36	 GFI20L14G42	 GFI20L14G50

* Not included in the INTEGRAL Guided surgery kits

FOR SLEEVES				
		Ø 3.6 mm ●	Ø 4.2 mm ●	Ø 5.0 mm ●
Guided step drills Ø 2.0/ 2.4 mm ●	Length 6.5 mm	 GFE24L6G36	 GFE24L6G42	 GFE24L6G50
	Length 8 mm	 GFE24L8G36	 GFE24L8G42	 GFE24L8G50
	Length 10 mm	 GFE24L10G36	 GFE24L10G42	 GFE24L10G50
	Length 12 mm	 GFE24L12G36	 GFE24L12G42	 GFE24L12G50
	Length 14 mm	 GFE24L14G36	 GFE24L14G42	 GFE24L14G50
Guided step drills Ø 2.4/ 3.0 mm ●	Length 6.5 mm	 GFE30L6G36	 GFE30L6G42	 GFE30L6G50
	Length 8 mm	 GFE30L8G36	 GFE30L8G42	 GFE30L8G50
	Length 10 mm	 GFE30L10G36	 GFE30L10G42	 GFE30L10G50
	Length 12 mm	 GFE30L12G36	 GFE30L12G42	 GFE30L12G50
	Length 14 mm	 GFE30L14G36	 GFE30L14G42	 GFE30L14G50
Guided step drills Ø 3.0/ 3.6 mm ●	Length 6.5 mm	-	 GFE36L6G42	 GFE36L6G50
	Length 8 mm	-	 GFE36L8G42	 GFE36L8G50
	Length 10 mm	-	 GFE36L10G42	 GFE36L10G50
	Length 12 mm	-	 GFE36L12G42	 GFE36L12G50
	Length 14 mm	-	 GFE36L14G42	 GFE36L14G50

FOR SLEEVES				
		Ø 3.6 mm	Ø 4.2 mm	Ø 5.0 mm
Guided step drills Ø 3.6/ 4.2 mm	Length 6.5 mm	-	-	 GFE42L6G50
	Length 8 mm	-	-	 GFE42L8G50
	Length 10 mm	-	-	 GFE42L10G50
	Length 12 mm	-	-	 GFE42L12G50
	Length 14 mm	-	-	 GFE42L14G50

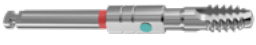FOR SLEEVES				
		Ø 3.6 mm	Ø 4.2 mm	Ø 5.0 mm
Guided taps Ø 3.4 mm	Short	 GTA34G36S	 GTA34G42S	 GTA34G50S
	Long	 GTA34G36L	 GTA34G42L	 GTA34G50L
Guided taps Ø 4.0 mm	Short	-	 GTA40G42S	 GTA40G50S
	Long	-	 GTA40G42L	 GTA40G50L
Guided taps Ø 4.6 mm	Short	-		 GTA46G50S
	Long	-		 GTA46G50L

FOR SLEEVES				
		Ø 3.6 mm	Ø 4.2 mm	Ø 5.0 mm
Guided implant mandrels BL	Short	 GMVG36S	 GMVG42S	 GMVG50S
	Long	 GMVG36L	 GMVG42L	 GMVG50L
Guided implant mandrels TL	Short	-	 GTIMG42S	 GTIMG50S
	Long	-	 GTIMG42L	 GTIMG50L
Guided implant wrenches BL	Short	 GCVG36S	 GCVG42S	 GCVG50S
	Long	 GCVG36L	 GCVG42L	 GCVG50L
Guided implant wrenches TL	Short	-	 GTIWG42S	 GTIWG50S
	Long	-	 GTIWG42L	 GTIWG50L

FIXATION OF THE SURGICAL GUIDE	REF.
Fixation pin	 OPFP20
Pin sleeves*	 OPFPSL20-3 (X3) OPFPSL20-5 (X5)
Pin drill	 OPFPD20

NON GUIDED INSTRUMENTS	REF.
Manual surgical wrench	 OPCS100
Reversible ratchet wrench	 INCC

* Not included in the INTEGRAL Guided surgery kits

FIRST DRILL GUIDED SURGERY

First drill guiding solution to ensure precise positioning of the implant axis and emergence point

- High precision instruments
- Compatibility with various kinds of implant planning software
- Compact kit and ergonomic instruments



FIRST DRILL SURGICAL KITS

FIRST DRILL Guided surgery kit



INMODIGM

FIRST DRILL Guided surgery kit - Empty

INMODIGMV

FIRST DRILL SETS

CONTENTS

FIRST DRILL Guided surgery set - Intro pack

1 x Guided centering drill
2 x Guided initial drills Ø 2.0 (S + L)

INSETIGM

FIRST DRILL Guided surgery set - Complete pack

1 x Guided centering drill
2 x Free guided pointer drills (S + XL)
3 x Guided initial drills M 2.0 (S + L + XL)

INSETIGMF

FIRST DRILL INSTRUMENTS

SHORT

LONG

XL

Guided pointer drill with stop Ø 2.0



OPGPDM

Guided pointer drill without stop Ø 2.0



OPGPDM S



OPGPDM XL

FIRST DRILL guided drill Ø 2.0



OPGDM-20S



OPGDM-20L



OPGDM-20XL

Pin drill



OPFPD20

ACCESSORIES

Guide sleeve Ø 2.0 mm



OPGS20
OPGS20-4 (X4)
OPGS20-10 (X10)

Pin sleeve



OPFPSL20-3 (X3)
OPFPSL20-5 (X5)

Fixation pin



OPFP20

Guided surgery sleeve insertion tool




OPGSTOOL-01

PROSTHETIC KIT

- All types of prothetic restorations
- Complete and compact
- Ergonomic: convenient to use
- Customizable: free gommets for additionnal instruments
- Secure: can be cleaned and sterilized



Axiom® Prosthetic kit

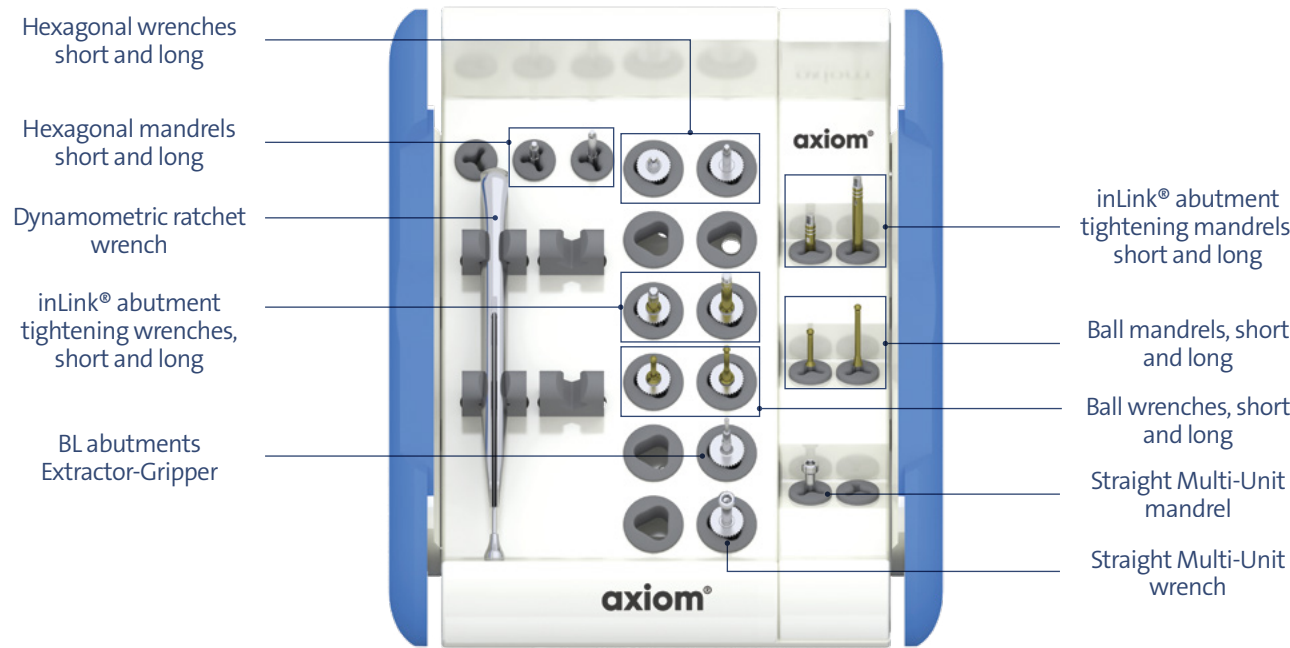

Prosthetic kit and instruments for Axiom® BL and Axiom® TL prosthesis

		USE
Axiom® Multi Level® prosthetic kit	 INMODOPP3	For Axiom® BL and TL prosthetic components
Empty Axiom® Multi Level® prosthetic kit	INMODOPP3V	
Update kit for Axiom® Multi Level® prosthetic kit	KITUOPP3	To update an Axiom® BL kit to an Axiom® ML kit
Dynamometrical ratchet wrench	 INCCD	15/25/35 N.cm

TORQ CONTROL®

Universal dynamometric torque wrench

Compatible with all implant systems

- Universal connection of the prosthetic mandrels
- 7 torque levels (10 to 35 N.cm), adapted to all protocols
- Automatic declutching when the selected torque is reached
- Micro-head, 100° angulated for easy access and optimal visibility, even in posterior sectors
- One-piece design for easy cleaning




More information p. 126

	USE	
Torq Control universal dynamometric wrench Including lubrication tip		Torques up to 35 N.cm Compatible with all implant systems
	15501*	

Instruments for Axiom® BL and Axiom® TL prosthesis

	SHORT	LONG	XL	USE
Hexagonal mandrel	 INMHECV 20 mm	 INMHELV 27 mm	 INMHEXLV* 35 mm	Main components on Axiom® BL and on Axiom® TL
Hexagonal wrench	 INCHECV 14 mm	 INCHELV 21.5 mm	 INCHEXLV* 29 mm	
Ball mandrel	 INBM100S INBM100S-3 (X3) 20.5 mm	 INBM100L INBM100L-3 (X3) 27 mm	 INBM100XL INBM100XL-3 (X3)* 35 mm	Angulated Access screw channels
Ball wrench	 INBW100S INBW100S-3 (X3) 23 mm	 INBW100L INBW100L-3 (X3) 29 mm	 INBW100XL INBW100XL-3 (X3)* 37 mm	
Angulated Access screw gripper	 AATool-M		 AATool AATool-4 (X4)	
Multi-Unit mandrel	 MUM100 16 mm	 MUM100L* 25.5 mm		For straight Multi-Unit abutments
Multi-Unit wrench	 MUW100			
Short prehensile tool for Multi-Unit	 MUWS*			
Wrench for analog for printed model 2 protective caps included PA-TOOL-01CAP	 PA-TOOL-01*			For all analogs for printed model
1st prosthetic set Mandrel: Hexagonal + Multi-Unit + Ball	PROTHSET-MS*	PROTHSET-ML*		
1st prosthetic set Wrench: Hexagonal + Multi-Unit + Ball	PROTHSET-WS*	PROTHSET-WL*		
Disconnecter and Gripper for Axiom® BL abutment	 INEXPS 32.4 mm	 INEXPL* 35.4 mm		
Harpoon abutment extractor for Axiom® BL abutments	 INEXPHS*	 INEXPHL*	 INEXPHK*	
Harpoon abutment extractor kit for Axiom® BL abutments	INEXPHSL-KIT containing INEXPHS + INEXPHL + INEXPHK			

*Not included in INMODOPP3 kit

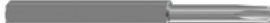
	SHORT	LONG	USE
Axiom® TL / mandrel	 TIM100S	 TIM100L	Axiom® TL
Axiom® TL / wrench	 TIW100S	 TIW100L	
AxiN® key Axiom® Multi Level®	 AXIN-TOOL-AML*		AxiN® for Axiom® BL and TL
Wrench for O'Ring	 INCOIO*		Dalbo ball abutments
Octogonal mandrel / O'Ring	 INMOICO*	 INMOILO*	
Dalbo® hand screw tightening spring	 O72609*		
LOCATOR® mandrel Ref. 8913	 O8913*		LOCATOR®
LOCATOR® key Ref. 8260	 O8260*		
LOCATOR® 3 in 1 key Ref. 8393	 O8393*		
Scan-Adapter tightening tool	 SATOOL-01*		All Scan Adapters

Laboratory kits

AXIN® LABORATORY KIT – AXIOM® MULTI LEVEL®		REF AXINLABKITAML		
1x AxIN® key Axiom® Multi Level®	 AXIN-TOOL-AML			
1x Ball wrench	 INBW100L			
6x AxIN® laboratory bases assorted (1 of each below references)				
Laboratory screw included	BONE LEVEL		TISSUE LEVEL	
GINGIVAL HEIGHT	Ø 4.0 mm	Ø 5.0 mm	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
1.5 mm	 AXIN152-27L41	 AXIN152-27L51	 AXIN156-02-L	 AXIN156-01-L
2.5 mm	 AXIN152-27L42	 AXIN152-27L52		

*Not included in INMODOPP3 kit

Digital Transfers for Intra Oral Scanners and analogs for printed models

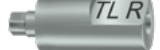SINGLE-UNIT RESTORATIONS ON AXIOM® BL AND AXIOM® TL			
	BONE LEVEL	TISSUE LEVEL	
Digital transfer – Single		PLATFORM R/4.8 mm	PLATFORM N/4.0 mm
Screw included	 152-27-DT 12.5 mm	 156-02-DT	 156-01-DT
Metal digital transfer - Single	 152-27-MDT 13.5 mm	 156-02-S-MDT	 156-01-S-MDT
Analog for printed model	 Laboratory screw included 152-27-PA 152-27-PA-4 (X4)	 156-02-PA 156-02-PA-4 (X4)	 156-01-PA 156-01-PA-4 (X4)
MULTIPLE-UNIT RESTORATIONS ON MULTI-UNIT ABUTMENTS			
	PLATFORM Ø 4.8 mm	PLATFORM Ø 4.0 mm	
Multi-Unit digital transfer			
Screw included	 151-03-DT-MU	 151-04-DT-MUN	
Multi-Unit metal digital transfer	 151-03-MDT*	 151-04-MDT*	
Analog for print model Multi-Unit	 151-03-PA 151-03-PA-4 (X4)	 151-04-PA 151-04-PA-4 (X4)	
Laboratory screw included			
MULTIPLE-UNIT RESTORATIONS ON AXIOM® TL			
	PLATFORM R/4.8 mm	PLATFORM N/4.0 mm	
Digital transfer (TL)			
Screw included	 156-02-DT-IL	 156-01-DT-IL	
Axiom® TL / analog for printed model	 156-02-PA 156-02-PA-4 (X4)	 156-01-PA 156-01-PA-4 (X4)	
Laboratory screw included			

Laboratory Scan-Adapters for Axiom® BL and TL Platforms


SINGLE-UNIT RESTORATIONS ON AXIOM® BL AND AXIOM® TL		
	BONE LEVEL	TISSUE LEVEL
Scan-Adapter		
Laboratory screw included	152-27-SAA	156-0X-SAA
MULTIPLE-UNIT RESTORATIONS ON MULTI-UNIT ABUTMENTS		
	MULTI-UNIT R/4.8 mm	MULTI-UNIT N/4.0 mm
Scan-Adapter		
Laboratory screw included	151-03-SAO	151-04-SAO
MULTIPLE-UNIT RESTORATIONS ON AXIOM® TL		
	TISSUE LEVEL	
Scan-Adapter		
Laboratory screw included	156-0X-SAO	
Scan-Adapter tightening tool		
	SA TOOL-01	

Scan-Adapters kits

AXIOM® MULTI LEVEL® SA KIT REF 151-02-KIT		
2x SAA Scan-Adapters Laboratory screw included	1 X AXIOM® BL  152-27-SAA	1 X AXIOM® TL  156-0X-SAA
4x Multi – Unit SAO Scan-Adapter Laboratory screw included	4 X MULTI-UNIT R/4.8 mm  151-03-SAO	
4x Axiom® TL SAO Scan-Adapter Laboratory screw included		4 X  156-0X-SAO
1 x Scan-Adapter screwdriver		1 X  SATOOL-01
AXIOM® SA UPGRADE KIT REF 151-02-KITUP		
1x Axiom TL single SAA Scan-Adapter Laboratory screw included		1 X  156-0X-SAA
4x Axiom® TL SAO Scan-Adapter Laboratory screw included		4 X  156-0X-SAO



Prosthetic screws for Axiom® BL

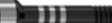		BONE LEVEL			
		TIGHTENING TOOL	BLACK	TITANIUM	LABORATORY
M1.6 prosthetic screw		 HEXAGONAL	 OPTS160 25 N.cm	 OPTS161 25 N.cm	Short  OPTS163 Manual Long  OPTS162 Manual
M1.6 screw for Axiom® Multi-Unit abutment			 OPMU160 25 N.cm		 OPMU161 OPMU161-4 (X4) Manual
X-Base® AA screws Axiom® BL	Gh 1.5 mm	 BALL	 OPFLEXS1-AA		 OPFLEXSL1-AA
	Gh 2.5 mm		 OPFLEXS2-AA		 OPFLEXSL2-AA
	Gh 3.5 mm		 OPFLEXS3-AA		 OPFLEXSL3-AA
AxIN® screw M1.6 Axiom® BL	H 1.5 mm		 AXIN152-27-S1 25 N.cm		 AXIN152-27-SL1 Manual
	H 2.5 mm		 AXIN152-27-S2 25 N.cm		 AXIN152-27-SL2 Manual

Prosthetic screws for Multi Unit on Axiom® BL and TL

	TIGHTENING TOOL	BLACK	TITANIUM	LABORATORY
M1.4 Multi-Unit screw	 HEXAGONAL	 MU140 15 N.cm	 MU140Z 15 N.cm	 MU141 Manual
				Short  MUT101 Manual Long  MUT102 Manual
Multi-Unit Angulated Access M1.4 prosthetic screw	 BALL	 MUAA140 15 N.cm	 MUAA141 15 N.cm	 MUAA142-4 (X4) Manual

Prosthetic screws for Axiom® TL

		TISSUE LEVEL				
		TIGHTENING TOOL	BLACK	TITANIUM	LABORATORY	
M1.6 Axiom® TL screw					Short 	Long 
Single-unit			TS160 25 N.cm	TS161 25 N.cm	TS163 Manual	TS162 Manual
Plural		HEXAGONAL				
			TS160P 25 N.cm	TS161P	TS163P-2 (X2)	TS162P-2 (X2)
X-Base® AA screws Axiom® TL	Single					
	Plural					
AxIN® screw M1.6 Axiom® TL						
Direct plural CAD/CAM restorations with AA on Axiom® TL						
			TSAA160P	TSAA161P	TSAA163P-2 (X2)	
Gripper for AA screws						
			AATOOL-M		AATOOL-4 (X4)	
						
			CR25-4		CR20-4	

AXIOM[®] 2.8

The solution for restricted mesiodistal spaces in the incisor area

- 3 implant lengths: 10, 12 and 14 mm
- 4 abutment angulations: 0°, 7°, 15° and 23°
- 4 gingival heights: 1.0, 2.5, 4.0 and 5.5 mm

A true narrow "two-part" implant

- Optimal management of treatment plan and soft tissues
- Same narrow 2.8mm diameter both for implant and abutment

COMPREHENSIVE PROSTHETIC RANGE

- PEEK temporary abutment
- One-stage or two-stage surgery



CONSTANT EMERGENCE PROFILE FROM HEALING PLUG TO FINAL ABUTMENT

- No undue tension on soft tissue



SMART IMPLANT PACKAGING

- Orange color code dedicated to Axiom[®] 2.8



DEDICATED SURGICAL & PROSTHETIC KIT

- Easy, practical and logic layout
- Safe Lock automatic prosthetic impactor



CONICAL CONNECTION


Tight seal morse tapered connection with 3° angle

A NARROW DIAMETER IMPLANT FOR THE INCISIVE SECTOR

For cases with narrow mesiodistal space, for single restorations.

STRAIGHT THREADED NECK

Designed for optimized primary bone anchorage

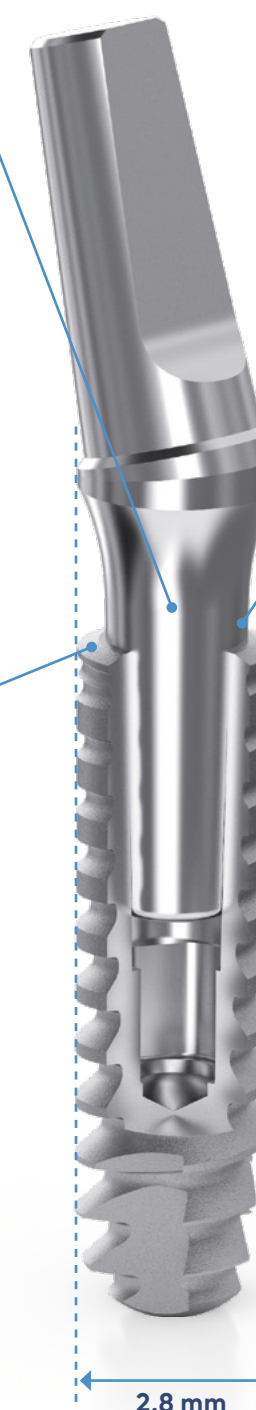

• AXIOM[®] 2.8

Cylindrical-conical shape implant with controlled and guided insertion.

Implants

STERILE Cover plug included OPIM028 LENGTH 10 mm 12 mm 14 mm	AXIOM [®] 2.8	
	Ø 2.8 mm	
		OP28100
		OP28120
		OP28140

Healing plugs in PEEK

HEIGHT	AXIOM [®] 2.8	
1.0 mm		OPHS210
2.5 mm		OPHS220
4.0 mm		OPHS240
5.5 mm		OPHS250

Temporary abutments in PEEK

HEIGHT	AXIOM [®] 2.8	
1.0 mm		OTP210
2.5 mm		OTP220
4.0 mm		OTP240
5.0 mm		OTP250

Impacted abutments

ANGULATION \ HEIGHT	AXIOM [®] 2.8			
	1.0 mm	2.5 mm	4.0 mm	5.5 mm
0°	 OPAT210	 OPAT220	 OPAT240	 OPAT250
7°	 OPAT211	 OPAT221	 OPAT241	 OPAT251
15°	 OPAT212	 OPAT222	 OPAT242	 OPAT252
23°	 OPAT213	 OPAT223	 OPAT243	 OPAT253

Impression and secondary components

	AXIOM [®] 2.8
Pop-in transfer	 OPPI028
Implant analog	 OPIA028
Abutment transfer	 OPTT028
Castable coping	 OPCA028
Protective cap	 OPPC028
Abutment analog	 OPAT028

Kit and instruments

	AXIOM [®] 2.8
Axiom [®] 2.8 kit	 INMODOP28
Safe Lock automatic prosthetic impactor Delivered with 3 impaction tips (1 x OPI100 + 1 X OPIP20 + 1 X OPIP400)	 6920
0° and 7° abutment tip	 OPIP100
15° and 23° abutment tip	 OPIP200
Definitive prosthetic impaction tip	 OPIP400
Round bur	 INFB20
Pointer drill	 OPPO15150
Initial drill Ø2.0 x 25.0	 OPFI20250
Drill Axiom [®] 2.8 Ø 2.6	 OPFE26250
Implant tap Ø 2.8	 OPTA28250

Instruments

	AXIOM [®] 2.8
Implant mandrel Axiom [®] 2.8	 OPMV028
Manual wrench Axiom [®] 2.8	 OPCV028
Depth gauge Axiom [®] 2.8 Ø 2.6	 OPJD026
Depth gauge Axiom [®] 2.8 Ø 2.0	 OPJD028
Threaded gripper wrench	 OPCF100
Prehensile wrench Axiom [®] 2.8	 OPOP028

Try-in abutments

		AXIOM [®] 2.8			
ANGULATION	HEIGHT	1.0 mm	2.5 mm	4.0 mm	5.5 mm
	0°	 OPAF210	 OPAF220	 OPAF240	 OPAF250
7°		 OPAF211	 OPAF221	 OPAF241	 OPAF251
15°		 OPAF212	 OPAF222	 OPAF242	 OPAF252
23°		 OPAF213	 OPAF223	 OPAF243	 OPAF253

MINI IMPLANT SYSTEM

One-piece mini implant with integrated Optiloc® prosthetic connection

- Usability: freedom of movement of optiloc® matrix system
- Flexibility: 3 matrix versions & 6 retention levels
- Ergonomics: smart dedicated surgical & prosthetic kits



INTEGRATED OPTILOC® PROSTHETIC CONNECTION

- Compact and slim connection
- Outstanding long lasting wear resistance due to the contact of the DLC* carbon-based coating of the Optiloc® connection with the retention insert PEEK material

TISSUE FAVORED DESIGN

- Soft tissue volume optimization

ONE-PIECE Ø2.6 MM MINI IMPLANT

- Minimally invasive treatment

CYLINDRICAL CONICAL SHAPE & SELF TAPPING THREAD

- Good primary stability

• ANTHOGRYR MINI IMPLANT SYSTEM

Indicated for patients with totally edentulous mandible and/or maxilla presenting severe reduced horizontal bone availability

Implants

STERILE	MINI IMPLANT
LENGTH	Ø2.6 mm Medical Grade V Titanium*
10 mm	 MIO26100
12 mm	 MIO26120
14 mm	 MIO26140

DRILLS	
Twist drill Ø1.5 mm Length: 17 mm Medical grade stainless steel	 INFH15170
Twist drill Ø2.0 mm Length: 17 mm Medical grade stainless steel	 INFH20170

SURGICAL INSTRUMENTS	
Gingival cutter Ø3.5 mm Medical grade stainless steel	 INDG035
Gauges Ø1.5 mm - Ø2.0 mm Depth and alignment gauges Medical Grade V Titanium*	 INJPAPF150200

MANDRELS AND WRENCHES	
Optiloc® mandrel Optiloc® adapter for contra-angle Medical grade stainless steel	 170.1-OPT - Reference to be ordered from Straumann
Optiloc® wrench Optiloc® adapter for ratchet wrench Medical grade stainless steel	 170.2-OPT - Reference to be ordered from Straumann
Ratchet wrench Medical grade stainless steel	 INCC
Surgical dynamometric wrench Medical grade stainless steel Reversible - Max. torque: 80 N.cm	 INCCDC


Instruments and surgical kit

DRILL STOPS				
Drill stops Ø2.0 mm Medical Grade V Titanium*				
	INB2004 Stopper for 4 mm drilling length	INB2010 Stopper for 10 mm drilling length	INB2012 Stopper for 12 mm drilling length	INB2014 Stopper for 14 mm drilling length
Drill stops kit Ø2.0 mm Medical Grade V Titanium* Kit including: Stopper for 4 mm drilling length (1 x INB2004) Stopper for 10 mm drilling length (1 x INB2010) Stopper for 12 mm drilling length (1 x INB2012) Stopper for 14 mm drilling length (1 x INB2014)	 KITINB			
MINI IMPLANT SURGICAL KIT				
Mini Implant surgical kit Complete kit including: 1 Gingival cutter (1 x INDG035) 1 Step drill Ø1.5 mm (1x INFH15170) 1 Step drill Ø2.0 mm (1 x INFH20170) 2 Gauges (1 x INJPAPF150200) 1 Ratchet wrench (1 x INCC) 1 Optiloc® adapter for contra-angle (1 x 170.1-OPT) 1 Optiloc® adapter for ratchet wrench (1 x 170.2-OPT)	 INMODMIOKIT			
Mini Implant surgical kit - Empty	INMODMIOV			

Optiloc® components

OPTILOC® IMPRESSION / MASTER CAST IMPRESSION	
Optiloc® forming / fixing matrix - White PEEK 4 pieces	 2102.0012-OPT - Reference to be ordered from Straumann
Optiloc® model analogue Aluminium 4 pieces	 2102.0024-OPT - Reference to be ordered from Straumann
OPTILOC® MATRIX HOUSINGS	
Optiloc® matrix housing Titanium 4 pieces	 2102.0001-OPT - Reference to be ordered from Straumann
Optiloc® matrix housing - Elliptic Titanium 4 pieces Recommended if extra-strong retention force of the matrix is required in the denture body	 2102.0009-OPT - Reference to be ordered from Straumann
Optiloc® matrix housing with attachment option Titanium 4 pieces <i>Special accessory</i> Used in situations requiring more retention. The attachment may be shortened according to the required height	 2102.0010-OPT - Reference to be ordered from Straumann
OPTILOC® RETENTION INSERTS	
Optiloc® retention insert - Red Retention force: extra-light Retention ≈ 300 g / ≈ 0.66 lbs* PEEK 4 pieces <i>Special accessory</i>	 2102.0003-OPT - Reference to be ordered from Straumann
Optiloc® retention insert - White Retention force: light Retention ≈ 750 g / ≈ 1.65 lbs* PEEK 4 pieces	 2102.0004-OPT - Reference to be ordered from Straumann
Optiloc® retention insert - Yellow Retention force: medium Retention ≈ 1 200 g / ≈ 2.65 lbs* PEEK 4 pieces	 2102.0005-OPT - Reference to be ordered from Straumann
Optiloc® retention insert - Green Retention force: strong Retention ≈ 1 650 g / ≈ 3.64 lbs* PEEK 4 pieces	 2102.0006-OPT - Reference to be ordered from Straumann
Optiloc® retention insert - Blue Retention force: extra-strong Retention ≈ 2 100 g / ≈ 4.63 lbs* PEEK 4 pieces	 2102.0007-OPT - Reference to be ordered from Straumann
Optiloc® retention insert - Black Retention force: ultra-strong Retention ≈ 2 550 g / ≈ 5.62 lbs* PEEK 4 pieces <i>Special accessory</i>	 2102.0008-OPT - Reference to be ordered from Straumann

* lbs = pounds

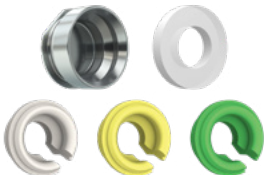OPTILOC® PROCESSING PACKAGE	
Optiloc® processing package Including: Optiloc® matrix housing - Titanium 2 pieces Optiloc® retention insert - White Light - 2 pieces Optiloc® retention insert - Yellow Medium - 2 pieces Optiloc® retention insert - Green Strong - 2 pieces Optiloc® mounting collar - 2 pieces Titanium / Silicone / PEEK	 5202.0001-OPT - Reference to be ordered from Straumann
OPTILOC® PROSTHETIC AIDS	
Optiloc® mounting collar Silicone 10 pieces	 2102.0011-OPT - Reference to be ordered from Straumann
Optiloc® processing spacer POM 4 pieces	 2102.0023-OPT - Reference to be ordered from Straumann
OPTILOC® TOOLS	
Optiloc® equipment box Including 3 tools: Optiloc® mounting & demounting tool for retention inserts - Brown 1 x 3202.0001-STM) Optiloc® mounting tool + model analogue reposition aid - Blue (1 x 3202.0002-STM) Optiloc® matrix housing extractor - Grey (1 x 3202.0003-STM)	 5102.0000-OPT - Reference to be ordered from Straumann
Optiloc® mounting & demounting tool for retention inserts - Brown Aluminium / Stainless steel	 3202.0001-OPT - Reference to be ordered from Straumann
Optiloc® model analogue reposition aid - Blue Aluminium / Stainless steel	 3202.0002-OPT - Reference to be ordered from Straumann
Optiloc® matrix housing extractor - Grey Aluminium / Stainless steel	 3202.0003-OPT - Reference to be ordered from Straumann

EQUIPMENT FOR IMPLANTOLOGY AND SURGERY


For more than 75 years, Anthogyr has put its know-how and expertise in the development and manufacturing of high precision surgical instrumentation to support dental surgeons and implantologists in their practice. Equipment, like implantology, is truly part of Anthogyr's history, with the same creative boldness to provide ingenious solutions with proven reliability.



TORQ CONTROL®

Universal torque wrench for controlled tightening

- 7 torque values adapted for all protocols
- High precision torque tightening whatever the implant connection
- Ergonomic design

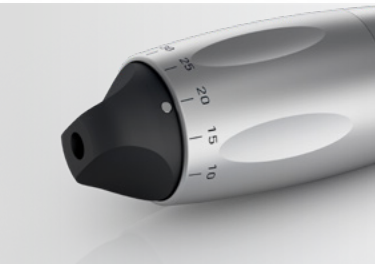

Compatible with all implant systems

- Universal connection
- 7 torque values (from 10 to 35 N.cm) adapted to all protocols



So safe

- High precision torque tightening
- Automatic declutching when the torque is reached
- Sterilizable in autoclave



Designed for efficiency

- Micro-head, 100° angulated
- For easy access and optimal visibility even in posterior sectors



Maximum comfort

- Only 130g for a better freedom of movement
- Manipulation 100% outside of mouth
- No lateral movement



Want to see more?



How to use Torq Control?



TORQ CONTROL®

Torq Control® universal torque wrench
Including lubrication tip



15501

XPERT UNIT®

Connected implantology motor

Ergonomically designed for high performance and efficiency.

- High torque capacity up to 80 N.cm
- Easy navigation and flexible programs
- Connected to a dedicated mobile app



XL touchscreen

- Ergonomics
- Working comfort and safety




Wireless pedal

- Freedom of movement
- Advanced foot navigation



Ultra-light micromotor and contra-angle

- Working comfort
- Dismantable contra-angle, compliant with the strictest maintenance protocols

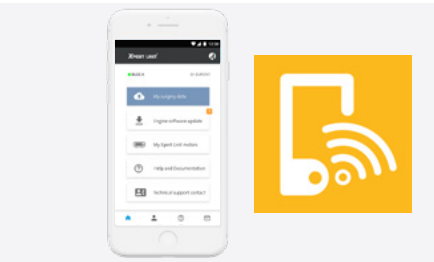

Connected to a dedicated mobile app

- New solutions for the traceability of surgery data
- Facilitated management of equipment
- And more friendly services...



Discover the app that makes your life easier

Please contact your sales representative to order the package bundle adapted to your needs

XPERT UNIT® PACK INCLUDING		
Unit + Wireless pedal + micromotor support + Fixation clips x 10 + bracket	XPERTUNIT	
MICROMOTORS		
Xpert Unit® micromotor without cable	1601242-001	1601244-001
Xpert Unit® micromotor cable	1601243-001	
CONTRA-ANGLES		
Short 20:1 Mont Blanc® implantology contra-angle	 16400XL	 16400X
Standard 20:1 Mont Blanc® implantology contra-angle	 10400XL	 10400X
Standard 20:1 Mont Blanc® implantology contra-angle with stop system	 10401XL	 10401X
ACCESSORIES		
Pedal charging adapter AC/DC 8W-5V power supply	XP216	(XP216UK for UK, Ireland, Malaysia and Singapore/XP216AU for Australia, New-Zealand)
Xpert Unit® power cable	XP217F - F type	(XP217G - G type for UK, Ireland, Malaysia and Singapore/ XP217I - I type for Australia, New-Zealand)
Xpert Unit® micromotor support	XP411	
Set of fixation clips for hanging irrigation line (x10)	XP412	
Xpert Unit® bracket	XP415	
CONSUMABLES		
Irrigation line for bag and bottle for Implanteo® and Xpert Unit® motors	32.F0209.00	
Screen protective film for Implanteo® and Xpert Unit® motors	32.U0022.00	
Motor protective sheath for Implanteo® and Xpert Unit® motors	22.U0001.00	

MONT BLANC® RANGE


Contra-angles and handpieces for surgery and implantology

- High performance
- Ergonomics
- Outcome of more than 75 years of experience



High torque resistance⁽¹⁾

- Capacity in all bone densities
- Contra-angle longevity

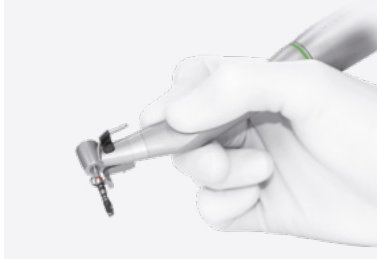

Elegance and ergonomics

- Streamlined and esthetic design
- Perfect angulation improving access to bone area
- Comfortable in the hand for precision drilling



Fully dismantlable

- Optimal decontamination
- Reliability and longevity



Stop system for high precision

- Adjustable stop system⁽²⁾ to determine the exact drilling depth
- Peace of mind

IMPLANTOLOGY CONTRA-ANGLES		
Mont Blanc® implantology contra-angle 20:1 SHORT*	 16400XL*	 16400X*
Mont Blanc® implantology contra-angle 20:1 STANDARD LENGTH	 10400XL	 10400X
Mont Blanc® implantology contra-angle 20:1 with stop system STANDARD LENGTH	 10401XL	 10401X

ACCESSORIES FOR MONT BLANC® 20:1 CONTRA-ANGLES		
Needle for internal spray for Mont Blanc® (supplied with 1 silicon tube)		10430
Sterilisable spray attachment (left side spray) for Mont-Blanc® implantology contra-angle		10475
Sterilisable spray attachment (right side spray) for Mont Blanc® implantology contra-angle		10475R
Single use spraying nozzles for Mont-Blanc® contra-angle (x10)		10476D
Complete stop holder + depth stop for Mont-Blanc®		10415
Mont Blanc® depth stop for ref. 10401X and 10401XL		10418
Clip for irrigation line for Mont-Blanc® contra-angle		2044
Contra-angle ruler		6458
Protection caps for optic fiber		10436
Y shaped stainless steel sterilizable coupling (delivered with 2 silicon tubes)		10420

SURGICAL HANDPIECES AND CONTRA-ANGLES		
Mont Blanc® surgical handpiece 1:1 straight	-	 12400X
Spraying nozzle for Mont Blanc® surgical handpieces		12430
Mont Blanc® surgical contra-angle 1:3	 12200XL	 12200X
Spraying nozzle for Mont Blanc® surgical contra-angles		12275

⁽¹⁾ Ref. 10400X, 10400XL, 16400X, 16400XL, 10401X, 10401XL
⁽²⁾ Ref. 10401X, 10401XL

*The short Mont Blanc® contra-angles are compatible with the Xpert Unit® and other motors with a nose shorter than 22.8 mm.

Axiom® implants packaging	128
Instructions for use	132
Choice of prosthetic component diameter	134
Try in abutments for Axiom® BL	135

USEFUL INFORMATION

IMPLANTS PACKAGING

Innovative, simple and reliable

IMPLANT « PROTECTED »
Titanium clamp

« NO TOUCH » SOLUTION
Direct pick-up of the implant with the contra-angle

ABILITY TO REPOSITION
the implant during surgery

INTUITIVE OPENING WITH ONE HAND

IMPLANT HANDLING USING « CLAMP »
Visible « PRESS » marking

CLOSURE SCREW
Packaging cap containing a cover screw for Axiom® BL and Axiom® TL implants and a cover plug for Axiom® 2.8

Implants packaging box

- Color coding «BL» and «TL»
- 4 sticky traceability labels
- Simple classification:
 - 1 color per diameter
 - 1 pictogram per implant
- Pre-cut box and tamper proof seal
- Sterile barrier

Color coding «TL» **PALE GREEN**

Color coding «BL» **PALE PINK**

Color coding						
Implant diameter	3.4 mm	4.0 mm	4.6 mm	5.2 mm	5.8 mm	6.4 mm

Color coding by implant diameter

en AXIOM IMPLANT Ø4.0 x 10.0
fr IMPLANT AXIOM Ø4.0 x 10.0
de AXIOM IMPLANTAT Ø4.0 x 10.0
it AXIOM IMPIANTO Ø4.0 x 10.0
es IMPLANTE AXIOM Ø4.0 x 10.0
pt IMPLANTE AXIOM Ø4.0 x 10.0

MD
CE
0459

REF OP40100
LOT XXXXX

YYYY-MM-DD
YYYY-MM-DD

ANTHOGYR ANTHOGYR ANTHOGYR
axiom BL REG axiom BL REG axiom BL REG
REF OP40100 REF OP40100 REF OP40100
LOT XXXXX LOT XXXXX LOT XXXXX

Supplied with 4 detachable labels

Packaging of prosthetic components

- Color coding «BL» and «TL»
- 2 color codes for three clear blisters
 - Blue: STERILE
 - Colorless: NON-STERILE
- Products supplied **STERILE**:
 - Fewer intra-oral handlings.
 - Intra-operative time saving.
- 4 detachable traceability labels*

Color coding «TL» **PALE GREEN**

Color coding «BL» **PALE PINK**

Color coding «Common to BL and TL, for ex MU» **WHITE**

Blue BLISTER (STERILE)
Examples:
standard titanium abutments, healing screws,
Multi-Unit abutments

Colorless BLISTER (NON-STERILE)
Examples: Flexibase®,
titanium temporary copings

Other colorless BLISTER (NON-STERILE)
Examples: transfers, analogs,
wrenches, drills.

* "Other colorless blister (non-sterile)" excluded.

Meaning of the symbols on blisters with traceability

Sealing lid

Do not re-use

Do not resterilize

Do not use if packaging damaged

CE marking - compliance with current regulations

Medical Device

Sterilized using irradiation

Consult instructions for use or consult electronic instructions for use

Manufacturer

Name and address of manufacturer

Date of manufacture

Use-by-date

Catalogue number

Batch code

ANTHOGYR

2237, Avenue André Lasquin
74700 SALLANCHES - FRANCE
MADE IN FRANCE

XXXX-MM-DD

XXXX-MM-DD

ifu.anthogyr.com

AXIOM IMPLANT Ø4.0 x 10.0

REF OP40100

LOT XXXXX

(01)XXXXXXXXXXXX(11)YYMMDD(17)YYMMDD(10)XXXXX

Example of repositionable traceability labels

- Easy, intuitive reading of the product references.
- Optimized component storage.

4 DETACHABLE LABELS

Product range

Reference

Batch number

Date of manufacture

Manufacturer

COMPONENT TYPE (DESIGNATION)

Article number

Batch number

Compliant with current regulations

Non sterile product

RECOMMENDED TIGHTENING TORQUE

UDI-DI code

COMPONENT PICTOGRAM

ASSOCIATED SCREW PICTOGRAM

ANTHOGYR

REF OPAT431

LOT XXXXX

Extractible-Removable

axiom BL

135°C

NON STERILE

Rx only

CE 0459

LLAFR03_00

de ASTHET TITAN-AUFBAU Ø4.0 15° H3

it MONCONE ESTET. TIT Ø4.0 15° H3

es PILAR EST. TIT. Ø4.0 15° ALT3

pt PILAR ESTÉ. TIT. Ø4.0 15° ALT3

(01)XXXXXXXXXXXX(11)YYMMDD(10)XXXXX

en AESTH. TIT. ABUT. Ø4.0 15° H3

fr MOIGNON ESTH. TIT. Ø4.0 15° H3

The abutment is extractible using the extractor INEXPS/L




EASY ACCESS TO INSTRUCTIONS FOR USE


You can access to all implants, prosthetic parts and instruments instructions for use on our on-line platform



Type the following address in your browser **ifu.anthogyr.com** or by clicking « **Instructions For Use** » on the Anthogyr website.

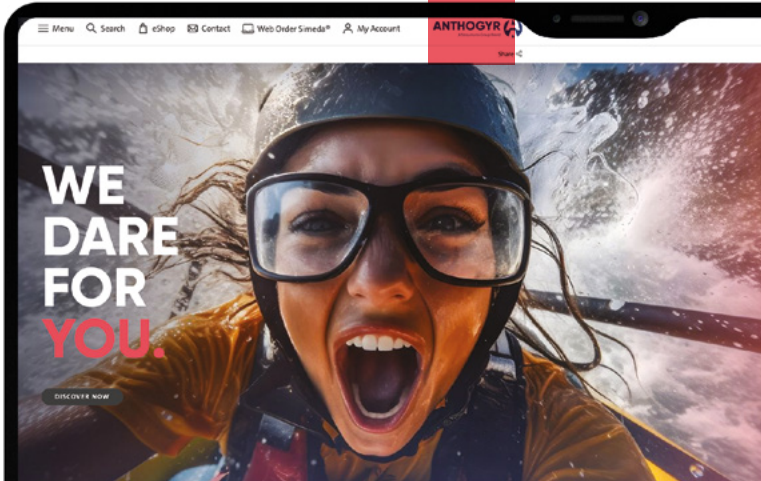


Select your country

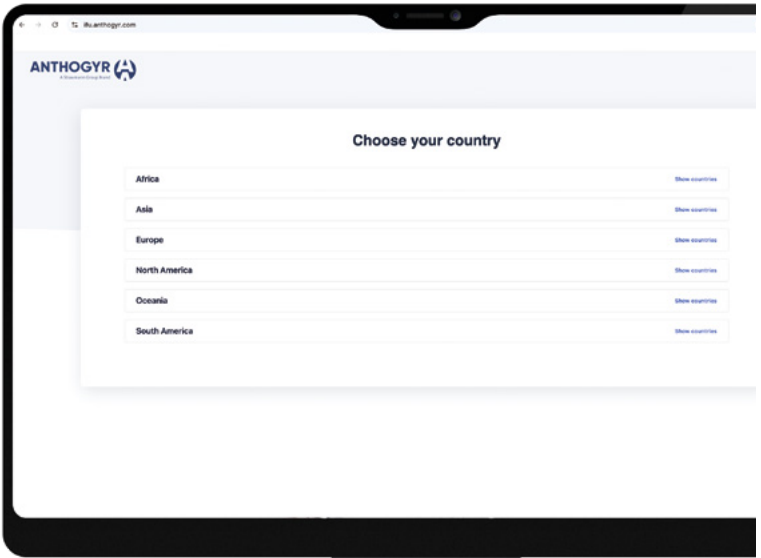

Enter the reference number, the product designation or the UDI-DI code and validate.

The IFU will open in a new window. Remember the IFU can be updated at any time. Always visit this website to get the latest version.

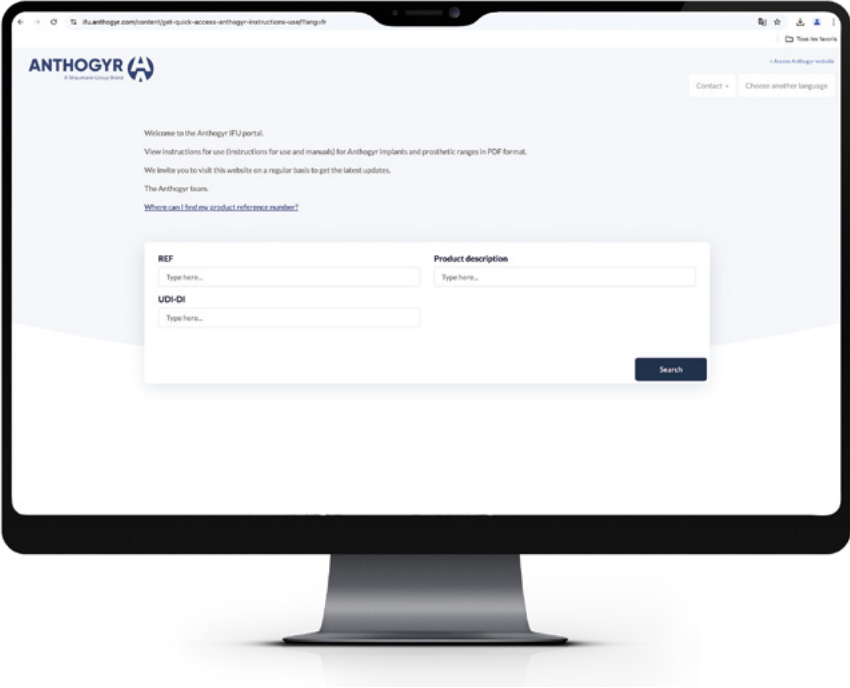

References, product designations and UDI-DI code are printed on packaging labels:

Product range

Reference

UDI-DI

ANTHOGYR REF OPST436 LOT XXXXX Extractible-Removable axiom BL	ANTHOGYR REF OPST436 LOT XXXXX Extractible-Removable axiom BL	ANTHOGYR REF OPST436 LOT XXXXX Extractible-Removable axiom BL
ANTHOGYR axiom BL MD 0459 STERILE R Rx only de TIT. STANDARD AUFBAU Ø4.0 H3/H6 it MONCONE STD TIT. Ø4.0 H3/H6 es PILAR TIT. ESTAND. Ø4.0 H3/H6 pt PILAR TITÂNIO STD Ø4.0 H3/H6 en STD.TITANIUM ABUT. Ø4.0 H3/H6 fr MOIGNON TITANE STD. Ø4.0 H3/H6		
(01)XXXXXXXXXXXX(11)YYMMDD(17)YYMMDD(10)XXXXXX		


CHOICE OF THE PROSTHETIC COMPONENT DIAMETER

Prosthetic diameter according to the tooth to be restored			RECOMMENDED EMERGENCE DIAMETER FOR AXIOM® BL			
MESIO-DISTAL WIDTH IN MM*			3.4	4.0	5.0	6.0
Maxilla	central incisor	7.6 - 10.5			●	●
	lateral incisor	5.3 - 8.3		●	●	
	canine	6.9 - 8.8		●	●	
	1 st premolar	6.0 - 8.2		●	●	
	2 ^d premolar	5.9 - 7.5		●	●	
	1 st molar	9.7 - 12.7			●	●
	2 ^d molar	8.7 - 11.4			●	●
Mandible	central incisor	4.7 - 6.2	●	●		
	lateral incisor	5.3 - 7.0	●	●		
	canine	6.0 - 8.1		●	●	
	1 st premolar	6.0 - 8.1		●	●	
	2 ^d premolar	6.4 - 8.8		●	●	
	1 st molar	9.7 - 12.5			●	●
	2 ^d molar	9.3 - 11.9			●	●

*Source: Lavergne, Bulletins et Mémoires de la Société d'Anthropologie de Paris, vol 1, serie XIII, 1974, 351-355.
● in first intention ● other possibility

TRY-IN ABUTMENTS FOR AXIOM® BL

Made of Medical titanium, with trilobe indexation

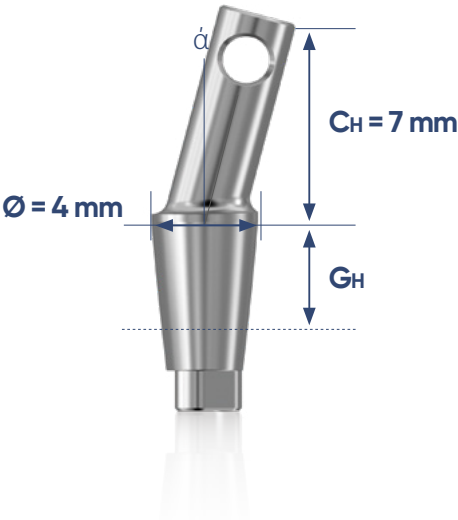

GINGIVAL HEIGHT	BONE LEVEL			
	STRAIGHT 0°	ANGULATION 7°	ANGULATION 15°	ANGULATION 23°
0.75 mm	 OPSF006	 OPAF00-7	 OPAF001	 OPAF002
1.5 mm	 OPSF016	 OPAF01-7	 OPAF011	 OPAF012
2.5 mm	 OPSF026	 OPAF02-7	 OPAF021	 OPAF022
3.5 mm	 OPSF036	 OPAF03-7	 OPAF031	 OPAF032
4.5 mm	 OPSF046	 OPAF04-7	 OPAF041	 OPAF042




Made in France

2025 Edition
C199_GB - 2026-07

Photo credits: Anthogyr, Getty images
All rights reserved - Product may differ from image shown.

Products availability may differ depending on the country.

www.anthogyr.com

FOLLOW US.
 **in**  **f**